



6200 series IT Carts

Model Numbers:

6200 series

Service and Parts Manual



FOR USE BY MIDMARK TRAINED TECHNICIANS ONLY

Table of Contents

General Info

GENERAL INFORMATION

<u>Symbols</u>	iii
<u>Ordering Parts</u>	iii
<u>Serial Number Location</u>	iii
<u>Specifications / Dimensions</u>	iv
<u>Password Retrieval (DC Cart software)</u>	v
<u>Lubricant for Rollerblade-style Casters</u>	v

Section A

TROUBLESHOOTING

not available

Section B

MAINTENANCE PROCEDURES

not available

Section C

ACCESS PROCEDURES

not available

Section D

WIRING DIAGRAMS

<u>DC Power Unit</u>	D-2
--------------------------------------	-----

Section E

EXPLODED VIEWS / PARTS LISTS

<u>6201 / 6202 / 6203 (Laptop / PC Carts)</u>	E-2
<u>6206 / 6207 / 6208</u>	
<u>(Kidney-style Laptop Carts)</u>	E-3
<u>6204 / 6205 (Small Laptop Carts)</u>	E-4
<u>6230 (Phlebotomy Carts)</u>	E-5
<u>6251 / 6252 / 6256 (DC Powered Carts)</u>	E-6

Digitally Linked Files

USER / SERVICE INFORMATION & FORMS

[Care & Cleaning Guide: all models](#)
[Cartridge Replacement: Non-Powered Carts](#)
[Installation / Troubleshooting Guide: DC Carts](#)
[Setup Guide: DC Carts](#)
[Operating Manual: DC Carts](#)
[User Guide: DC Carts](#)

<u>Parts Order Form</u>	004-0755-00
<u>Comments Form</u>	004-0756-00
<u>Color Selector</u>	www.midmark.com

(*) Indicates multiple pages due to model / serial number break(s).

Symbols



DANGER

Indicates an imminently hazardous situation which will result in serious or fatal injury if not avoided. This symbol is used only in the most extreme conditions.



WARNING

Indicates a potentially hazardous situation which could result in serious injury if not avoided.



Caution

Indicates a potentially hazardous situation which may result in minor or moderate injury if not avoided. It may also be used to alert against unsafe practices



Equipment Alert

Indicates a potentially hazardous situation which could result in equipment damage if not avoided.

Note

Used for special instructions or additional information.

Ordering Parts

The following information is required when ordering parts:

- Serial number & model number
- Part number for desired part
(Refer to Section E: Exploded Views & Parts Lists)

Non-warranty parts orders may be faxed to Midmark using the Fax Order Form in the back of this manual.

For warranty parts orders, call Midmark's Technical Service Department with the required information.


Hours: 8:00 am to 5:00 p.m. EST (Monday thru Friday)

Phone: 1-800-Midmark (1-800-643-6275)

Serial Number Location

Non-Powered Models

(6201 / 6202 / 6204 / 6205 / 6206 / 6207 / 6230)

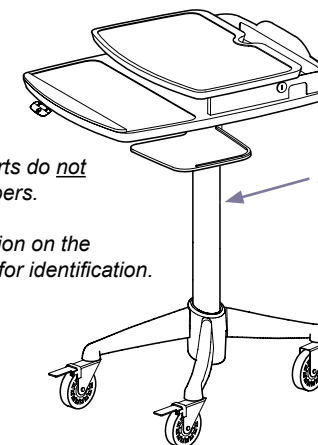


Non-powered carts do not have serial numbers.

Use the information on the Product ID label for identification.

Product ID Label

DC Powered Models

(6251 / 6252 / 6256)

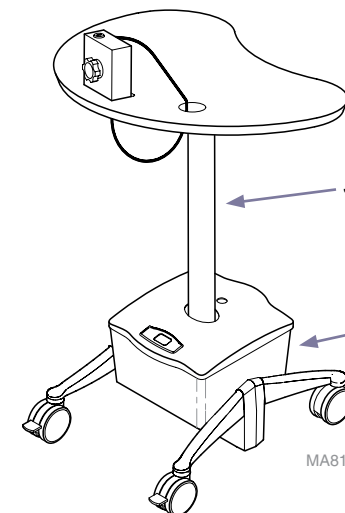

Serial Number



MA81611

Specifications / Dimensions

Weight Capacity:	
models 6201 / 6202 / 6203 / 6206 / 6208 / 6230 6251 / 6252 / 6256	30 lbs (13.5 kg)
models 6204 / 6205 / 6207	10 lbs (4.5 kg)
Base Dimensions:	
Non-Powered Workstations: models 6201 / 6202 / 6204 / 6205 / 6206 / 6207	18" x 20" (46 cm x 51 cm)
Phlebotomy Workstations: model 6230	18" x 20" (46 cm x 51 cm)
DC Powered Workstations: models 6251 / 6252 / 6256	18" x 22" (46 cm x 56 cm)
Worksurface Dimensions:	
models 6201 / 6202 / 6251 / 6252	20" D x 20" W (51 cm x 51 cm)
models 6206 / 6207 / 6256	22" D x 32" W (56 cm x 81 cm)
models 6230	26" D x 20" W (66 cm x 51 cm)
models 6204	15" D x 26" W (38 cm x 66 cm)
models 6205	18" D x 28" W (46 cm x 71 cm)
Regulatory Compliance:	
Non-Powered Workstations: models 6201 / 6202 / 6204 / 6205 / 6206 / 6207	ANSI / BIFMA
Phlebotomy Workstations: model 6230	ANSI / BIFMA
DC Powered Workstations: models 6251 / 6252 / 6256	ANSI / BIFMA UL 60601-1 CAN/CSA-C22.2 No. 601-1-M90

Laptop Workstations

(models: 6201 / 6202 / 6204 / 6205 / 6206 / 6207)

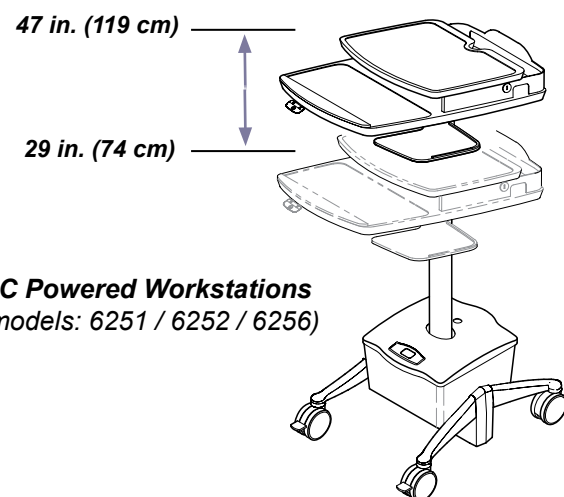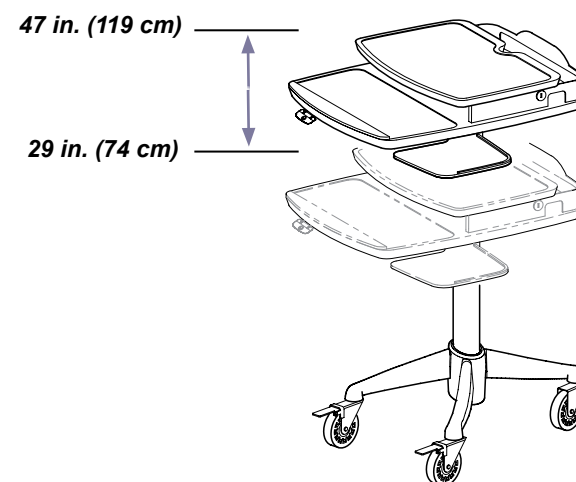


DC Powered Workstations

(models: 6251 / 6252 / 6256)

Phlebotomy Workstation

(model: 6230)

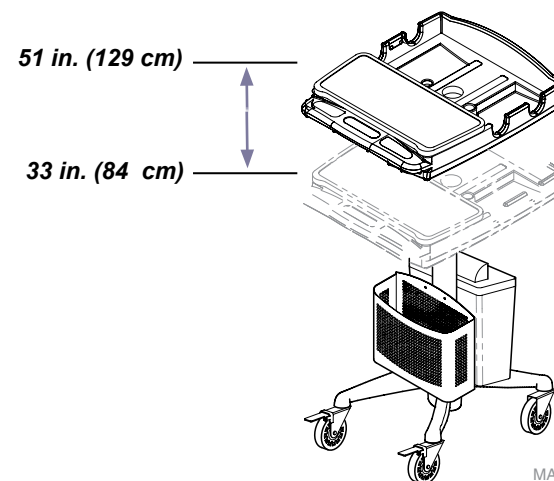

MA8158i

Password Retrieval (DC Cart Software)


The factory default password is “SERIES 301”. This password will always work, and gives access to all advanced functions.

If password is forgotten:

1. From the Nurses screen, click the “Advanced” button.
2. Type in “SERIES 301” (hit Enter).
3. Click “Setup” from menu bar, then select “GUI Setup”.
4. Click the “Setup Password” button.
5. Enter a new password, then click “Save”
6. Click the “Exit” button.
7. Click “File” on the menu bar, then select “Software Exit”.
(*The GUI program will close completely.*)
8. Restart the GUI program by clicking the windows “Start” button.
9. Select “Programs” / “Startup”, then click on “Series 301”.
10. Once the Nurses screen is operating, the new password for Advanced functions should be active

Lubricant for Rollerblade-style Casters

To lubricate a noisy caster, use standard bearing grease or WD40.

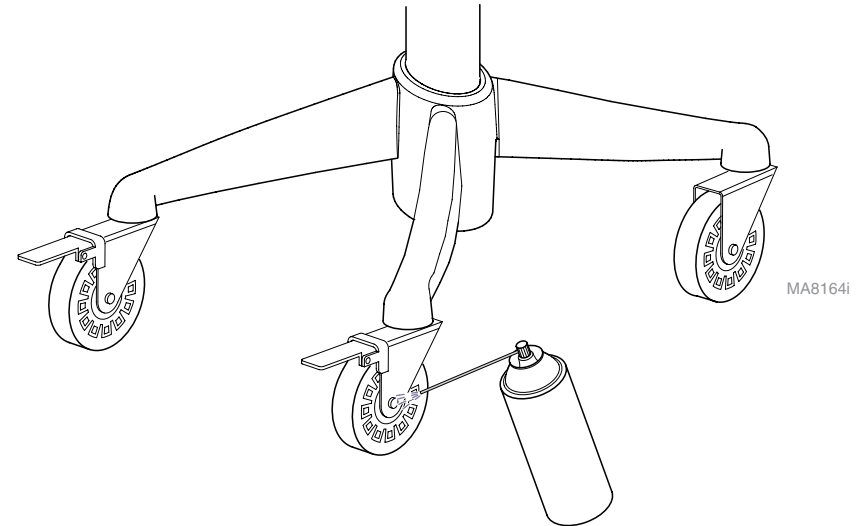



Models:

Serial Numbers:

Wiring Diagrams

Section D

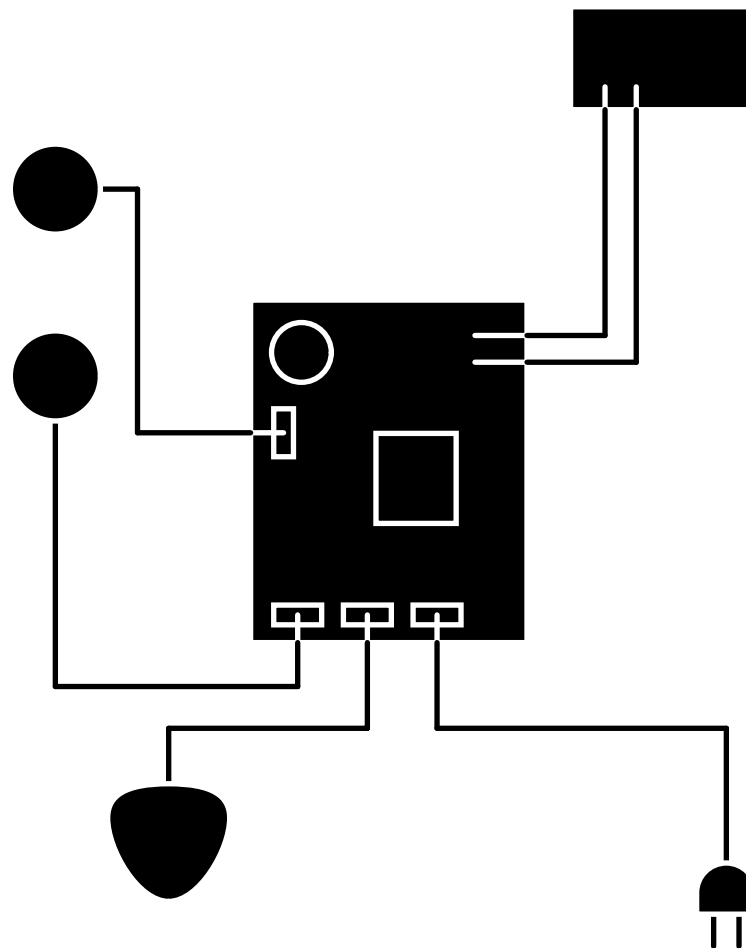

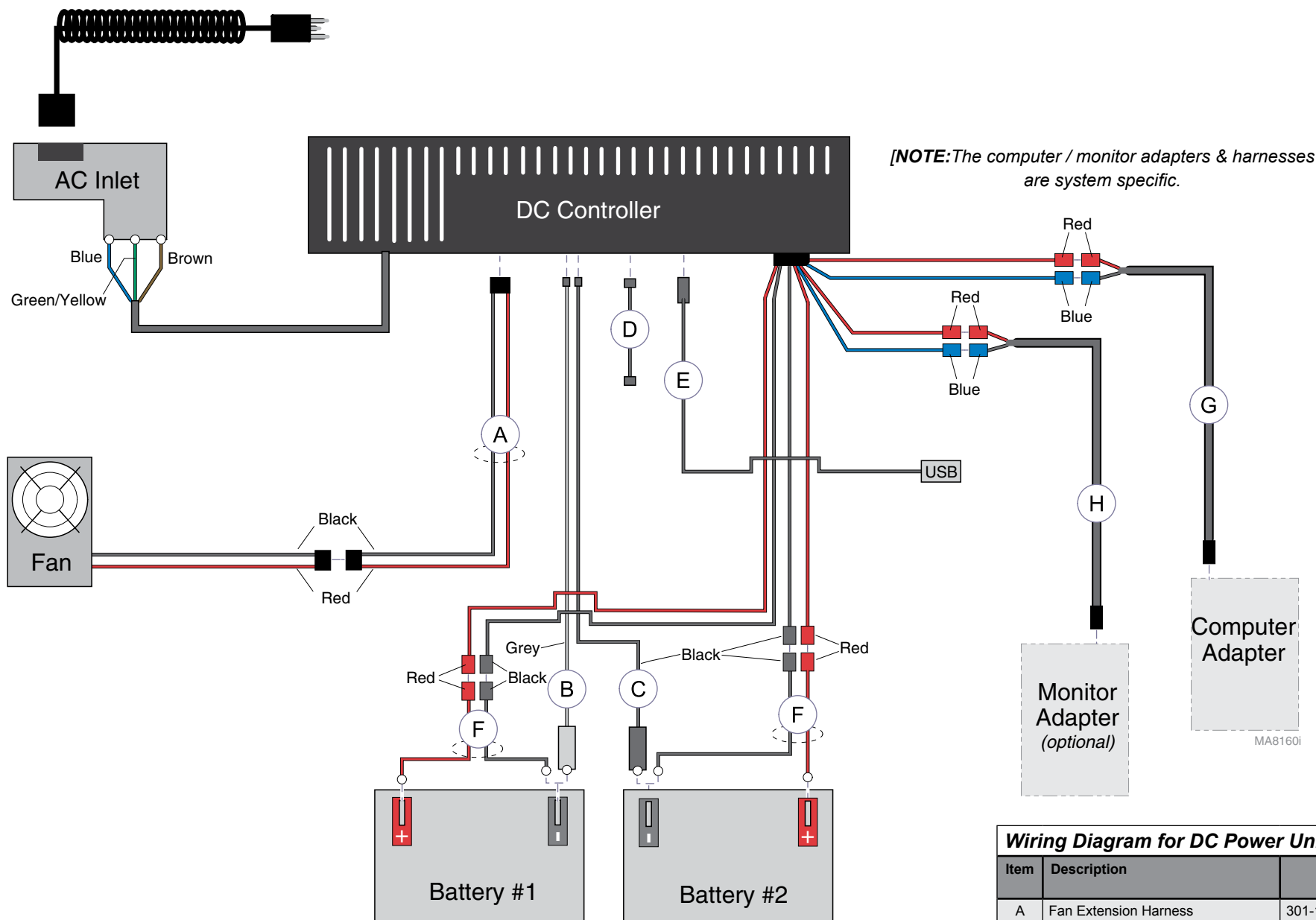

Wiring Diagram for DC Power Unit

Item	Description	Part Number
A	Fan Extension Harness	301-1005-WH-004
B	Battery #1 Temp. Sensor Harness	304-1001-AS-001
C	Battery #2 Temp. Sensor Harness	304-1001-AS-002
D	CAT 5 Cable	ACC-0110
E	USB Harness	CBL-0046
F	Red/Black Battery Harness (qty: 2)	305-1001-WH-008
G	Computer Adapter Harness	-
H	Monitor Adapter Harness (optional)	-

Always verify model & serial number


Models: 6251 / 6252 / 6256

Serial Numbers: all

Section E

Exploded Views & Parts Lists

Laptop Carts	E-2
Kidney-style Laptop Carts	E-3
Phlebotomy Carts	E-4
DC Powered Carts	E-5

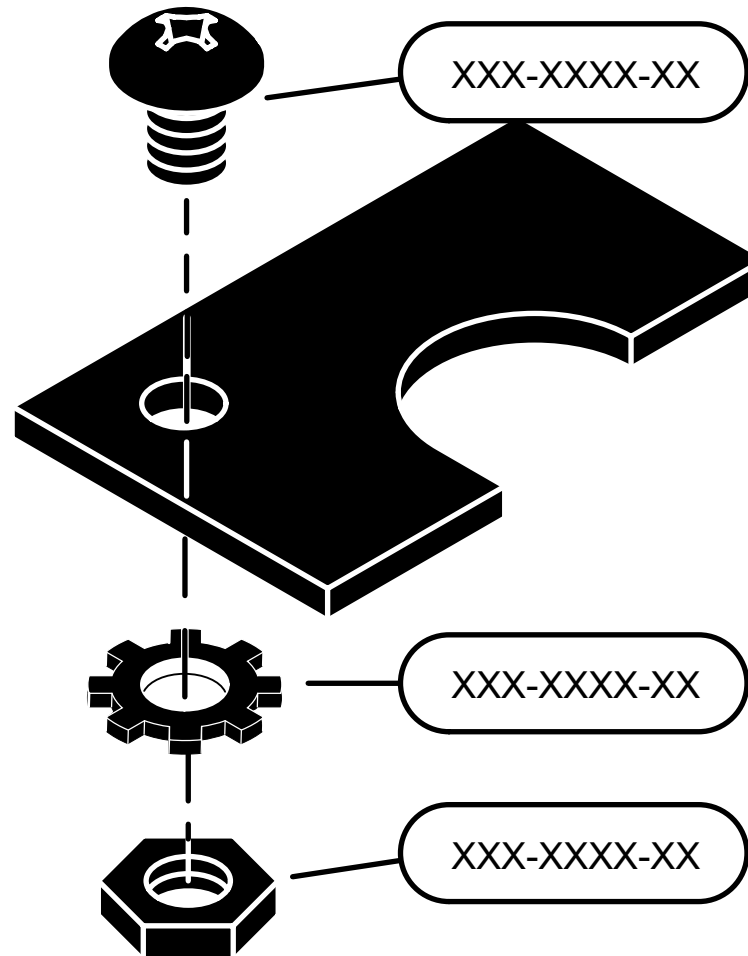

Laptop / PC Carts

[Flat Panel Monitor Mount](#) E-33


[Flat Worksurface](#) E-10

[Laptop / Monitor Mount Worksurface](#) E-14


[Swing Arm Mouse Pad](#) E-24

[Base / Cartridge Components \(Non-Powered Carts\)](#).....E-7


6202

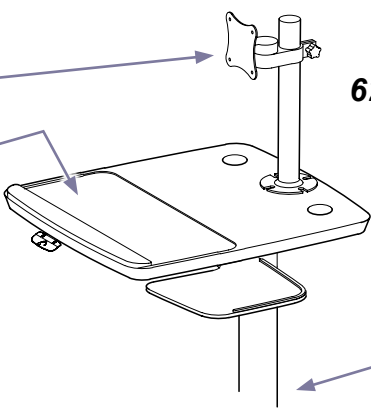


6201

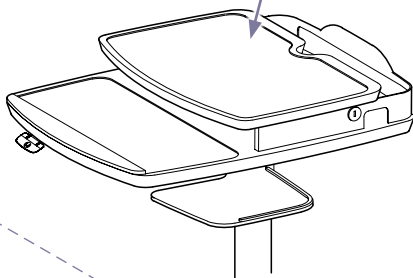

6203

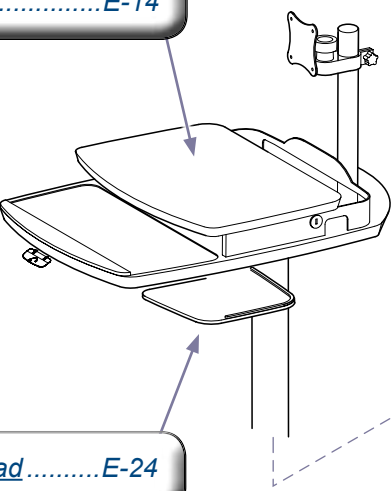

MA8135i

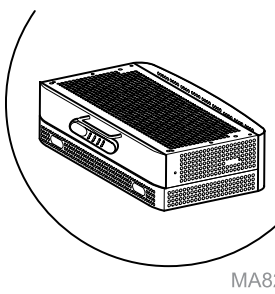

MA8247i

[Drawer / Bin Options](#) E-36


[Slim Line CPU Holder](#).....E-25

[Laptop Worksurface](#) E-13

- Accessories:**
- [Coiled Power Cord \(9A435\)](#)..... E-17
 - [Retractable Cable Reel \(9A436\)](#)..... E-18
 - [Cable Mgt Holder \(9A437\)](#)..... E-19
 - [Utility Basket \(9A439\)](#)..... E-21
 - [Multi-Purpose Shelf \(9A440\)](#) E-22
 - [Bar Code Holder \(9A441\)](#)..... E-23
 - [CPU Holder \(9A444\)](#)..... E-26
 - [Power Distribution Strip \(9A445\)](#) E-27
 - [IV Pole \(9A448\)](#) E-28
 - [Sharps Container \(9A450\)](#)..... E-29
 - [Lockable Drawer \(9A451\)](#)..... E-30
 - [Cable Mgt Holder \(9A452\)](#)..... E-31
 - [Cable Mgt Holder \(9A453\)](#)..... E-32
 - [Rear Basket \(9A467\)](#)..... E-35

Models:	6201 / 6202 / 6203	
Serial Numbers:	all	

Kidney-style Laptop Carts

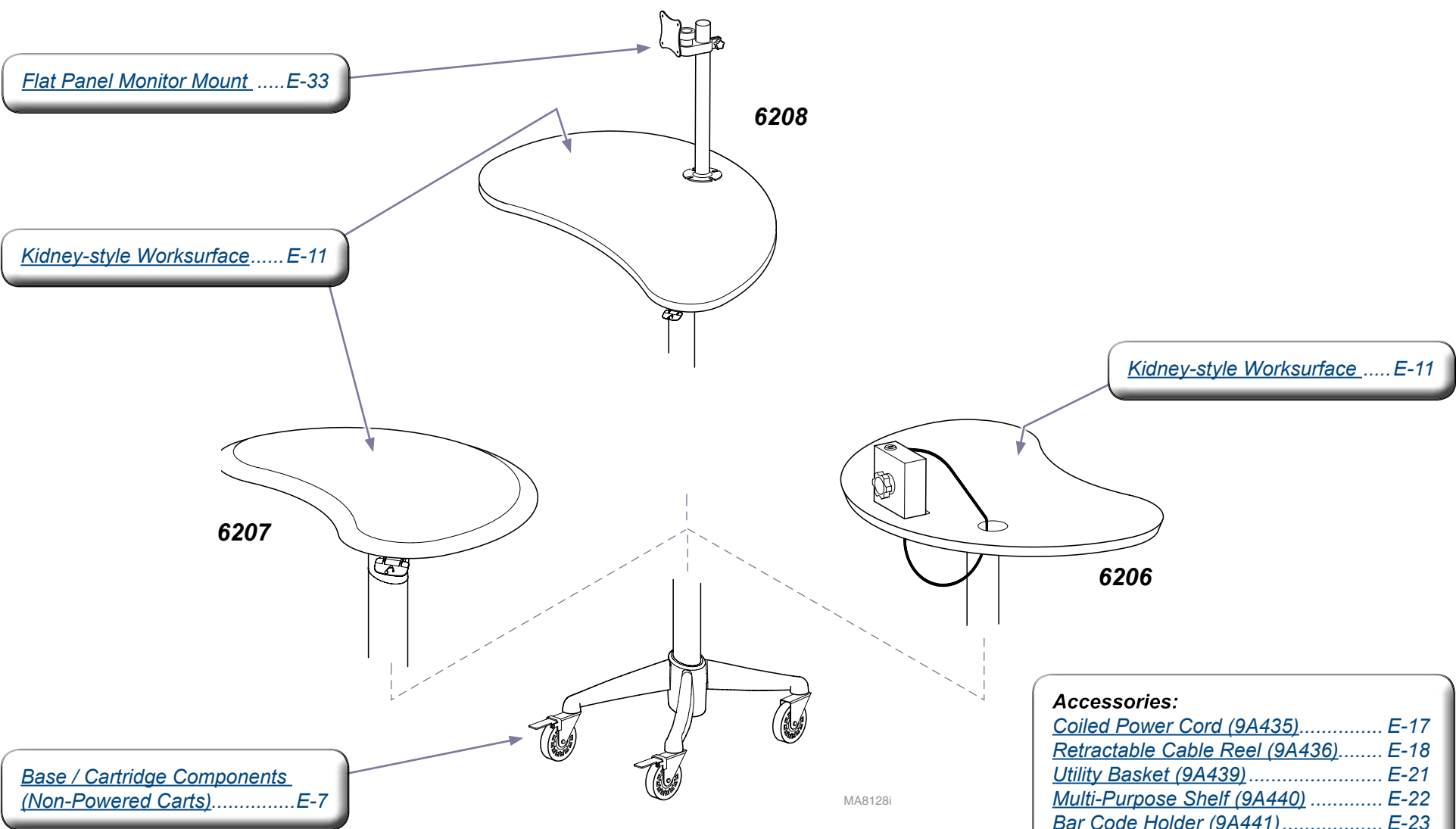

Accessories:


Coiled Power Cord (9A435)	E-17
Retractable Cable Reel (9A436)	E-18
Utility Basket (9A439)	E-21
Multi-Purpose Shelf (9A440)	E-22
Bar Code Holder (9A441)	E-23
Slim Line CPU Holder (9A443)	E-25
CPU Holder (9A444)	E-26
Power Distribution Strip (9A445)	E-27
IV Pole (9A448)	E-28
Sharps Container (9A450)	E-29
Cable Mgt Holder (9A452)	E-31
Cable Mgt Holder (9A453)	E-32

Models:	6206 / 6207 / 6208	
Serial Numbers:	<i>all</i>	

Small Laptop Carts

[Small Kidney-style
Worksurface E-12](#)

[Small Mayo-style
Worksurface E-12](#)

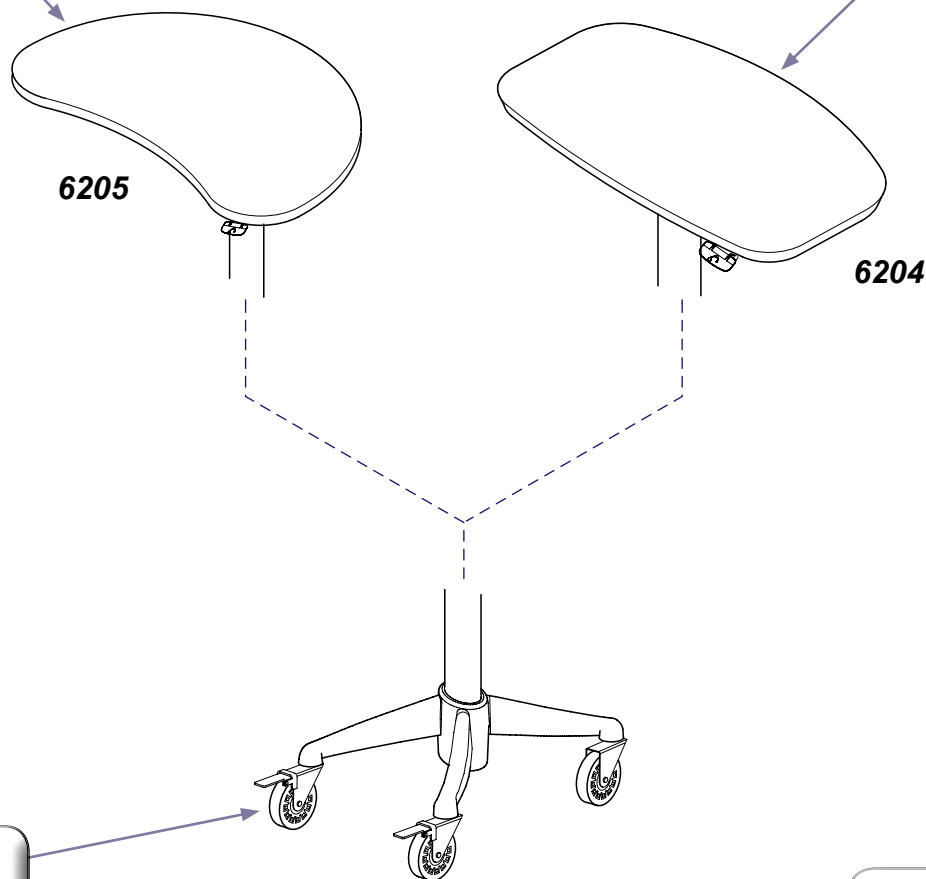

[Base / Cartridge Components
\(Non-Powered Carts\).....E-7](#)

MA8245i

Accessories:

[Retractable Cable Reel \(9A436\)..... E-18](#)
[Utility Basket \(9A439\)..... E-21](#)
[Slim Line CPU Holder \(9A443\)..... E-25](#)
[CPU Holder \(9A444\)..... E-26](#)
[IV Pole \(9A448\)..... E-28](#)
[Sharps Container \(9A450\)..... E-29](#)
[Utility Basket \(9A466\)..... E-34](#)


Models:	6204 / 6205	
Serial Numbers:	<i>all</i>	

Phlebotomy Carts

Accessories:
[Keyless Security System \(9A438\) E-20](#)

[Phlebotomy Worksurface Components E-15](#)

[Utility Basket Accessory E-21](#)

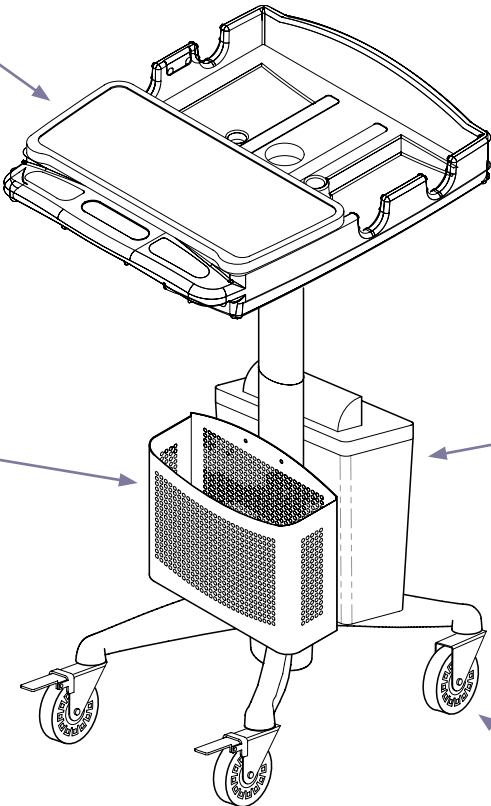

[Sharps Container Accessory E-29](#)

[Base / Cartridge Components \(Non-Powered Carts\) E-7](#)

MA8129i

Models:	6230	
Serial Numbers:	<i>all</i>	

DC Powered Carts

[Laptop Worksurface Components](#) E-13

[Flat Panel Monitor Mount Accessory](#) E-33

[Flat Worksurface Components](#) E-10

[DC Power Unit](#) E-9

[Drawer / Bin Options](#) E-36

[Swing Arm Mouse Pad Accessory](#) E-24


[Kidney-style Worksurface Components](#) E-11

[Base / Cartridge Components \(DC Carts\)](#) E-8


- Accessories:**
- [Cable Mgt. Holder \(9A437\)](#) E-19
 - [Utility Basket \(9A439\)](#) E-21
 - [Multi-Purpose Shelf \(9A440\)](#) E-22
 - [Slim Line CPU Holder \(9A443\)](#) E-25
 - [IV Pole \(9A448\)](#) E-28
 - [Sharps Container \(9A450\)](#) E-29
 - [Cable Mgt Holder \(9A452\)](#) E-31
 - [Cable Mgt Holder \(9A453\)](#) E-32
 - [Rear Basket \(9A467\)](#) E-35

Models:	6251 / 6252 / 6256	
Serial Numbers:	<i>all</i>	

Optional Caster Leg Assemblies

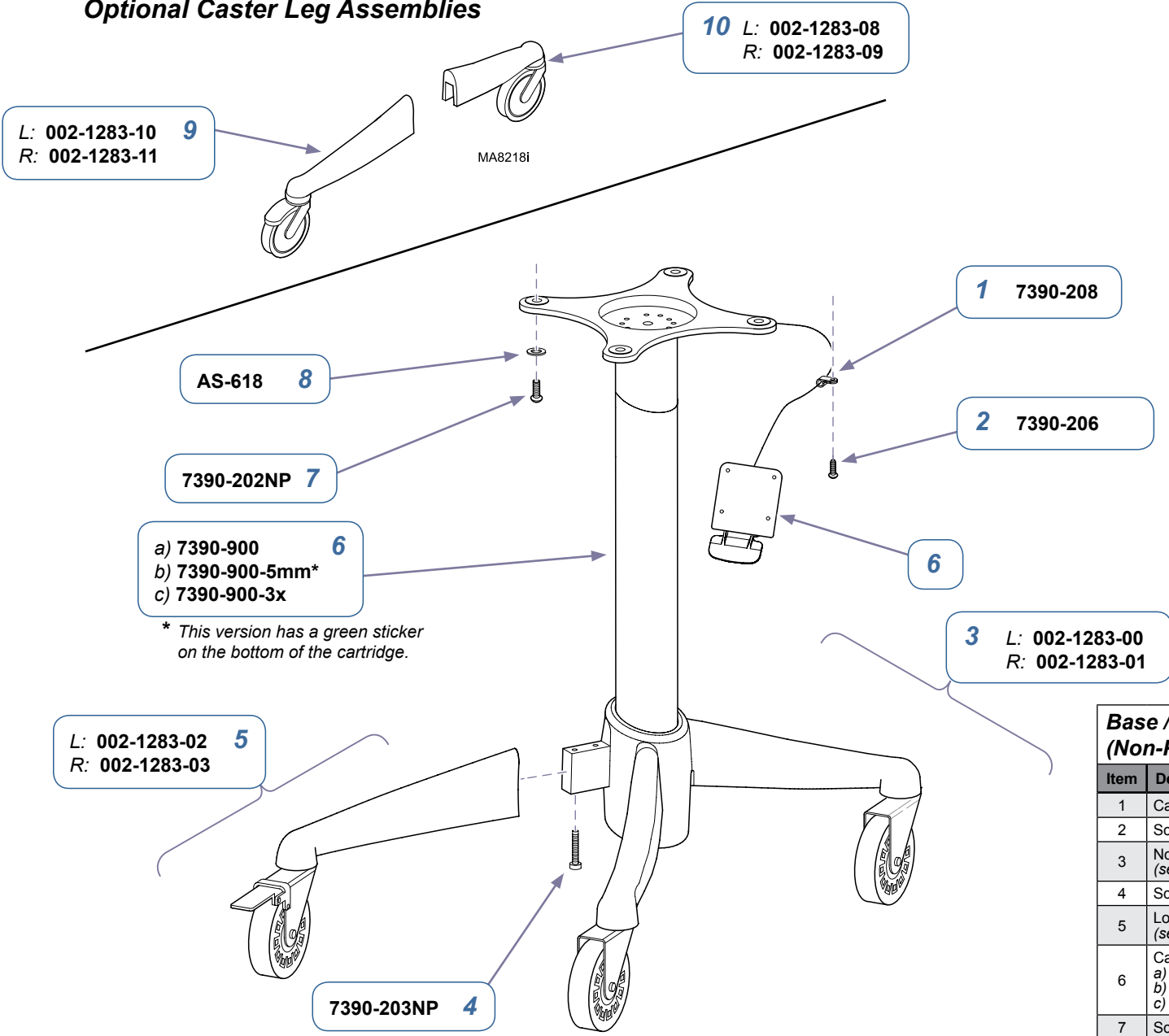

Caster / Leg Assemblies
(bottom view)

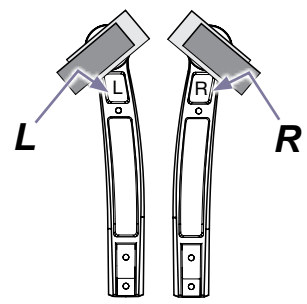

Base / Cartridge Components (Non-Powered Carts)		
Item	Description	Qty.
1	Cable Clamp	1
2	Screw (#8 x 1/2")	5
3	Non-Locking Caster / Leg Assembly (see bottom view detail illustration)	2
4	Screw (1/4-20 x 1 3/4")	8
5	Locking Caster / Leg Assembly (see bottom view detail illustration)	2
6	Cartridge Assembly (includes cable release lever) a) for models: 6201 / 6202 / 6203 / 6206 / 6230 b) for model: 6204 / 6205 / 6207 c) for model: 6203-001 <u>only</u>	1
7	Screw (3/8-16 x 7/8")	4
8	Lockwasher	4
9	Optional Locking Caster Leg Assembly (see bottom view detail illustration)	2
10	Optional Non-Locking Caster Leg Assembly (see bottom view detail illustration)	2

Models:	6201 / 6202 / 6203 / 6204 / 6205 / 6206 / 6207 / 6230
Serial Numbers:	all

AS-618 12

7390-202NP 11

1 models 6251 & 6256: 7390-810
model 6252: 7390-810-3X

F6NO.75BK 13

2 AS-487

3 AS-821

4 AS-488

5 L: 002-1283-08
R: 002-1283-09

6 L: 002-1283-04
R: 002-1283-05

7 7390-203NP

8 L: 002-1283-06
R: 002-1283-07

9 L: 002-1283-10
R: 002-1283-11

10

MA8144i

Caster / Leg Assemblies (bottom view)

Base / Cartridge Components (DC Powered Carts)


Item	Description	Qty.
1	Cartridge Assembly (includes cable release lever)	1
2	Column Clamp	2
3	Stafast JCM Nut	4
4	Threaded Rod	2
5	Current Production Units Non-Locking Caster Leg Assembly (see bottom view detail illustration)	
6	Earlier Production Units Non-Locking Caster Leg Assembly (see bottom view detail illustration)	2
7	Screw (1/4-20 x 1 3/4")	4
8	Earlier Production Units Locking Caster Leg Assembly (see bottom view detail illustration)	
9	Current Production Units Locking Caster Leg Assembly (see bottom view detail illustration)	2
10	Refer to: " DC Power Unit "	Ref
11	Screw (3/8-16 x 7/8")	4
12	Lockwasher	4
13	Sleeving	1

Models: 6251 / 6252 / 6256

Serial Numbers: all

E-8


Always verify model & serial number

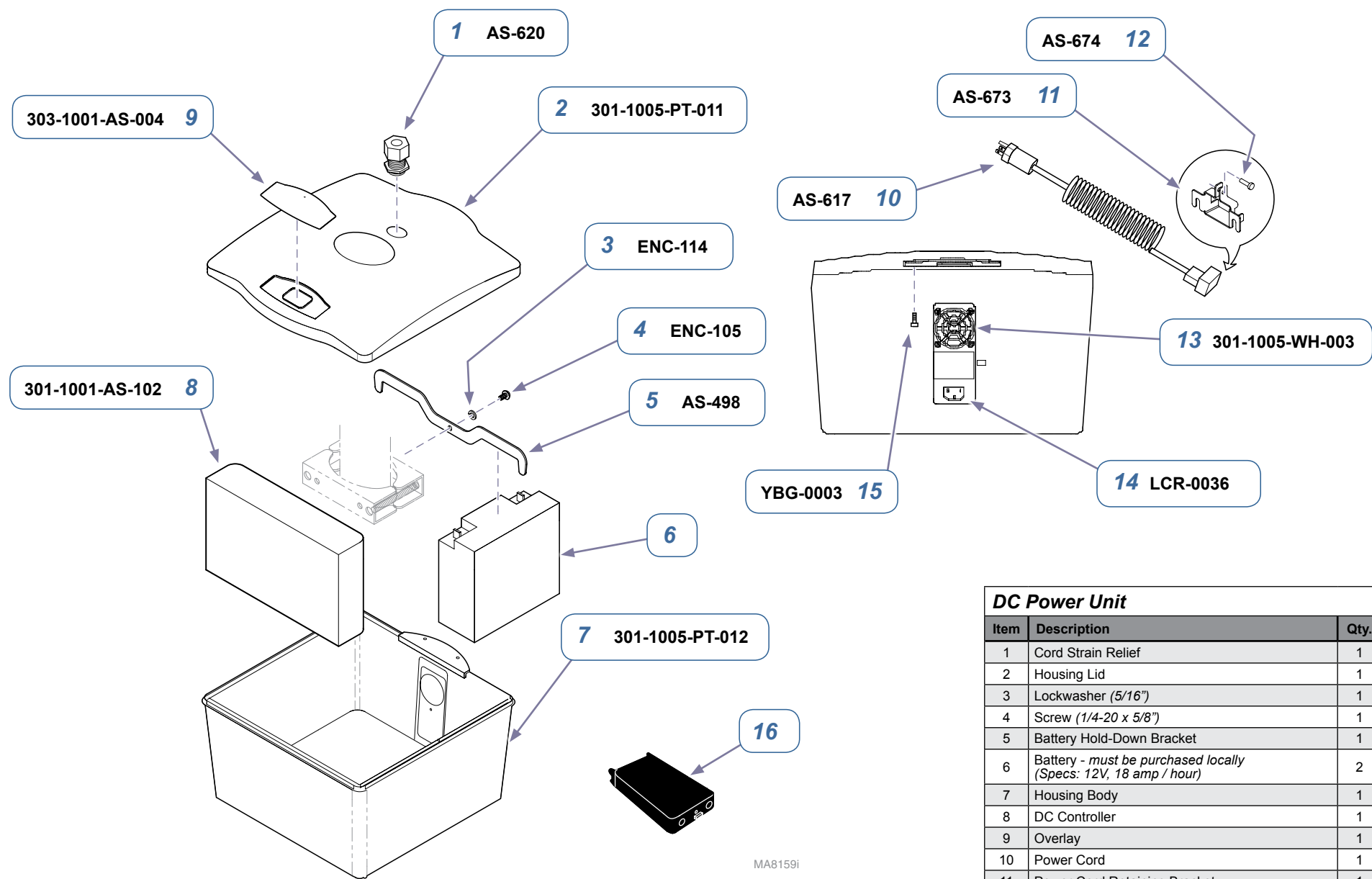

DC Power Unit		
Item	Description	Qty.
1	Cord Strain Relief	1
2	Housing Lid	1
3	Lockwasher (5/16")	1
4	Screw (1/4-20 x 5/8")	1
5	Battery Hold-Down Bracket	1
6	Battery - <i>must be purchased locally</i> (Specs: 12V, 18 amp / hour)	2
7	Housing Body	1
8	DC Controller	1
9	Overlay	1
10	Power Cord	1
11	Power Cord Retaining Bracket	1
12	Screw (#4-40 x 3/8")	1
13	Fan	1
14	AC In-Line Filter	1
15	Screw (8-32 x 3/8")	3
16	Brick - <i>Must contact Lind Electronics for support and replacement. Call 1-800-697-3702 or E-mail techsupport@lindelectronics.com</i>	1

Models:	6251 / 6252 / 6256	
Serial Numbers:	<i>all</i>	


Color Selector..... www.midmark.com

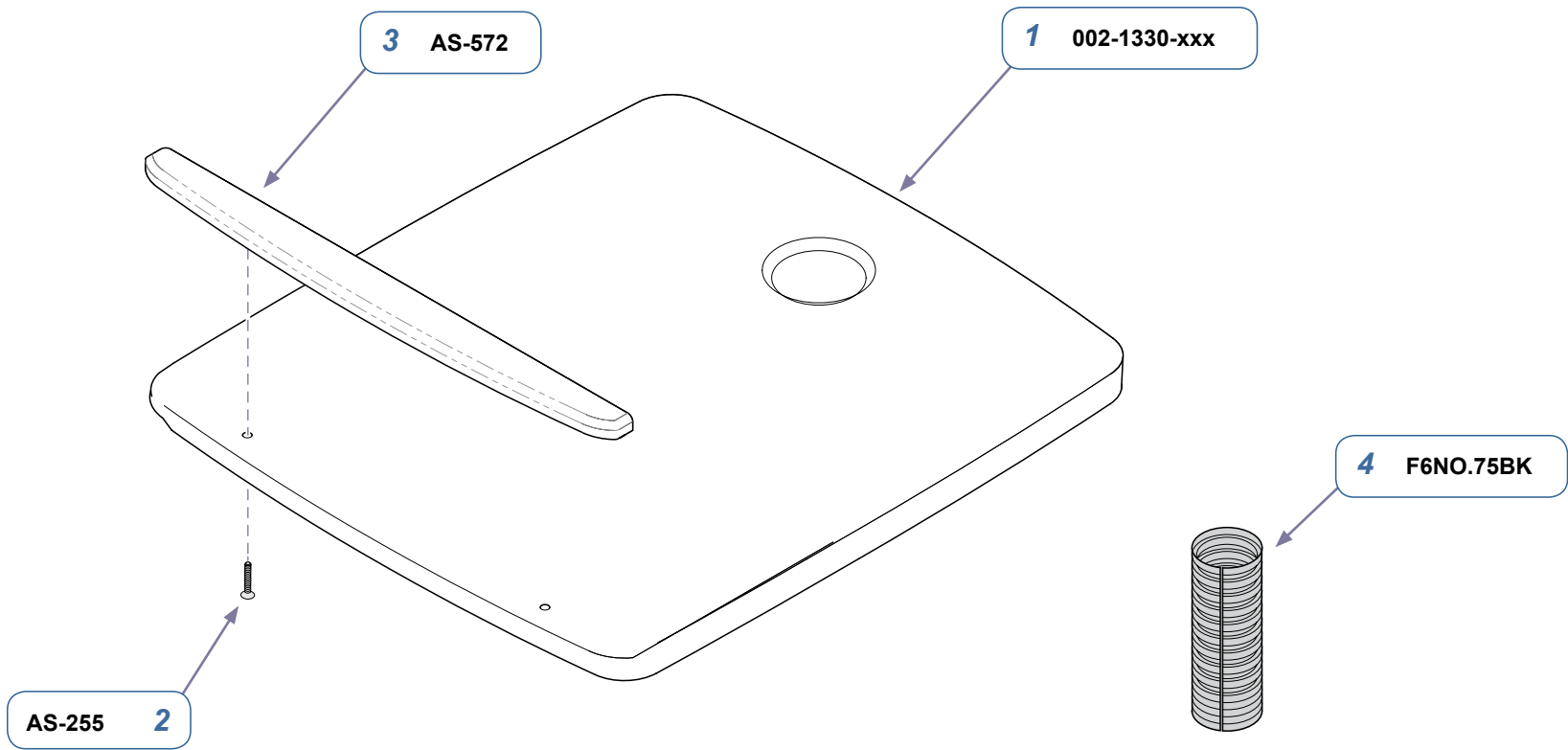

MA8134i


Models:	6202 / 6252	
Serial Numbers:	<i>all</i>	

Flat Worksurface Components		
Item	Description	Qty.
1	Laptop Worksurface Assembly (<i>incl. items 2 & 3</i>)	1
2	• Screw (#8-32 x 1")	2
3	• Wrist Rest	1
4	Sleeving	1

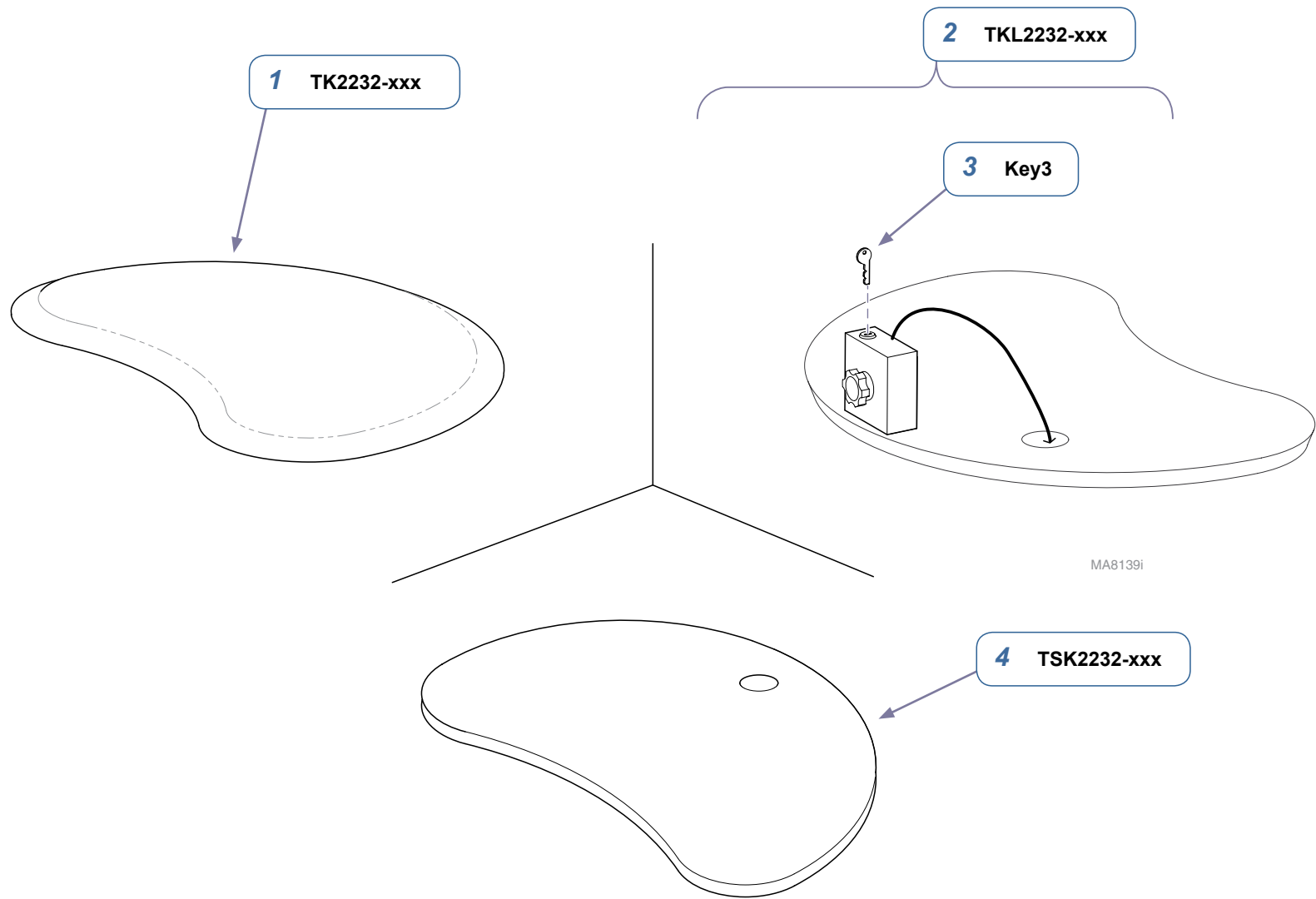

Models:	6206 / 6207 / 6208 / 6256	
Serial Numbers:	all	


Kidney-style Worksurfaces		
Item	Description	Qty.
1	Kidney-style Worksurface	1
2	Kidney-style Worksurface w/ Cable Lock (includes item 3)	1
3	• Key	1
4	Kidney-style Worksurface w/ monitor mount hole	1

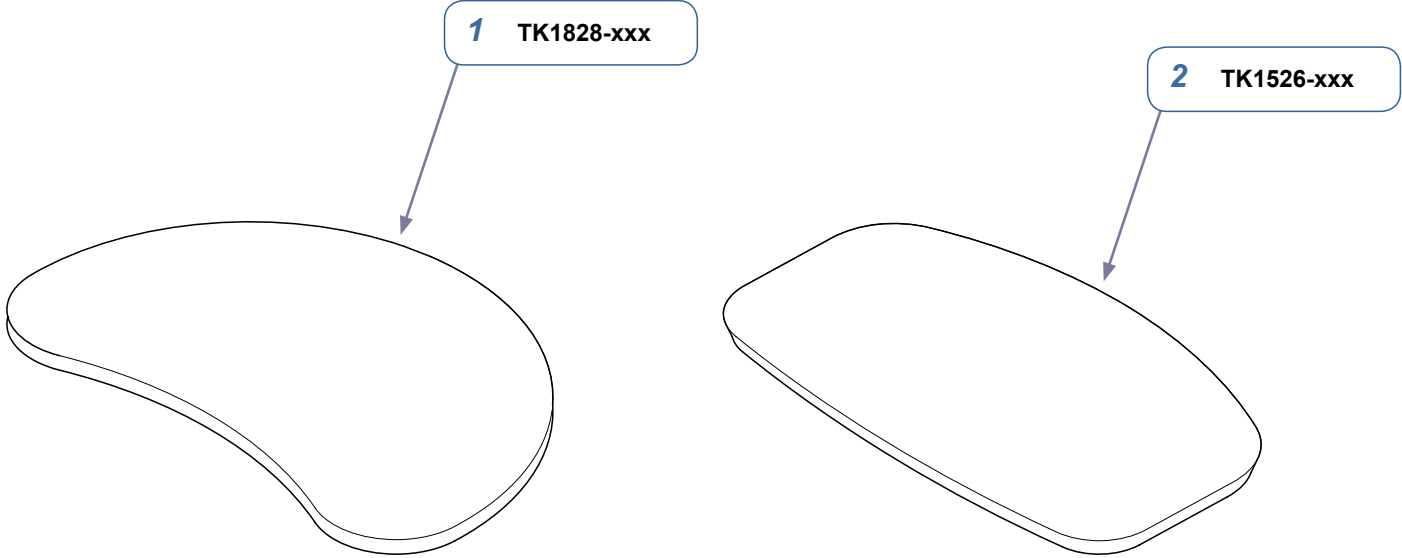

MA8246i

Models:	6204 / 6205	
Serial Numbers:	<i>all</i>	


Small Laptop Worksurfaces		
Item	Description	Qty.
1	Small Kidney-style Worksurface	1
2	Small Mayo-style Worksurface	1

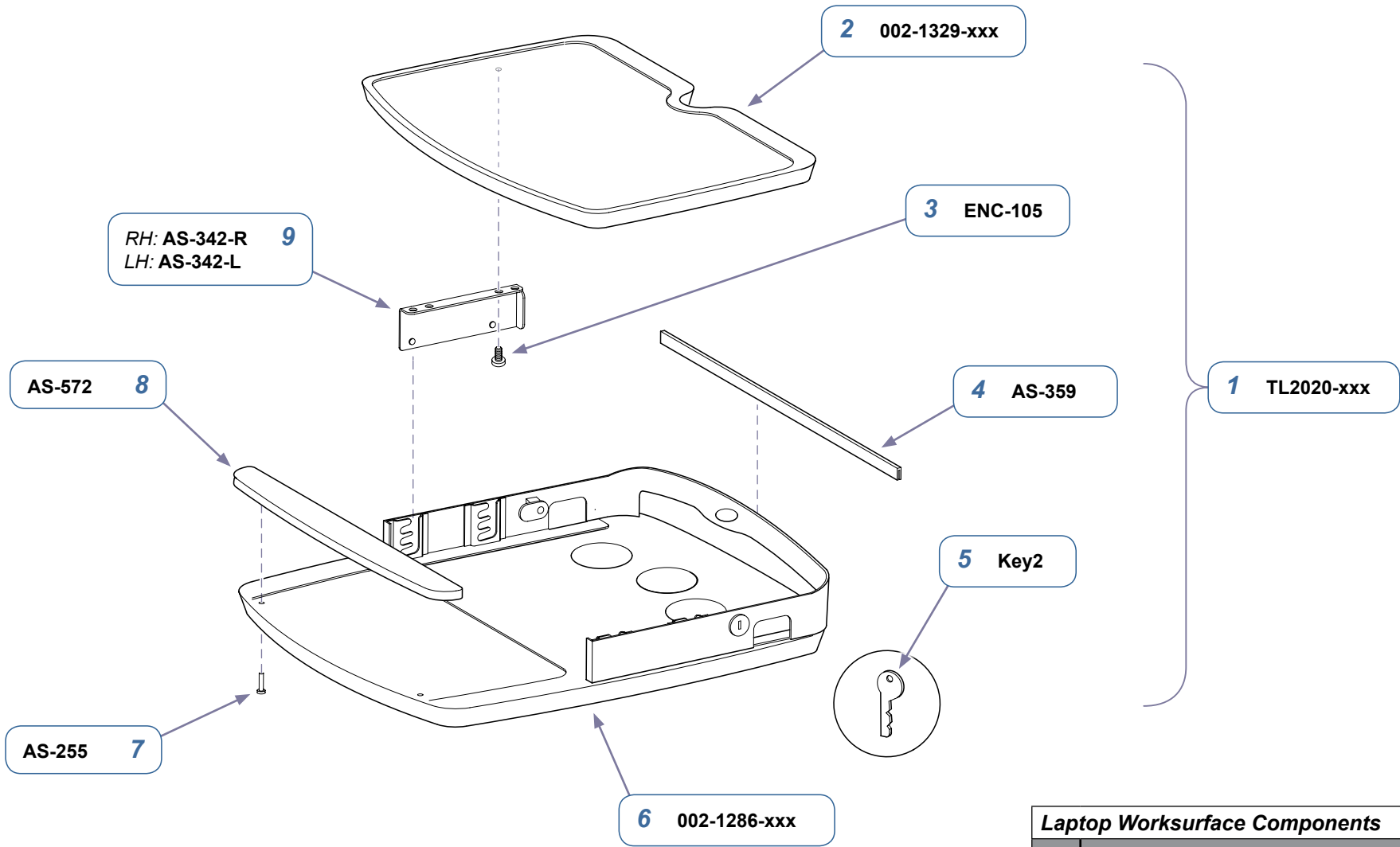

MA8140i


Laptop Worksurface Components		
Item	Description	Qty.
1	Laptop Worksurface Assembly (incl. 2 thru 9)	1
2	• Upper Worksurface	1
3	• Screw (1/4-20 x 5/8")	4
4	• Vinyl Edge Trim (per foot)	2 ft
5	• Key	1
6	• Lower Worksurface / Lock Assembly	1
7	• Screw (#8-32 x 1")	2
8	• Wrist Rest	1
9	• Upper Mounting Bracket (LH shown)	2

Models:	6201 / 6251	
Serial Numbers:	all	

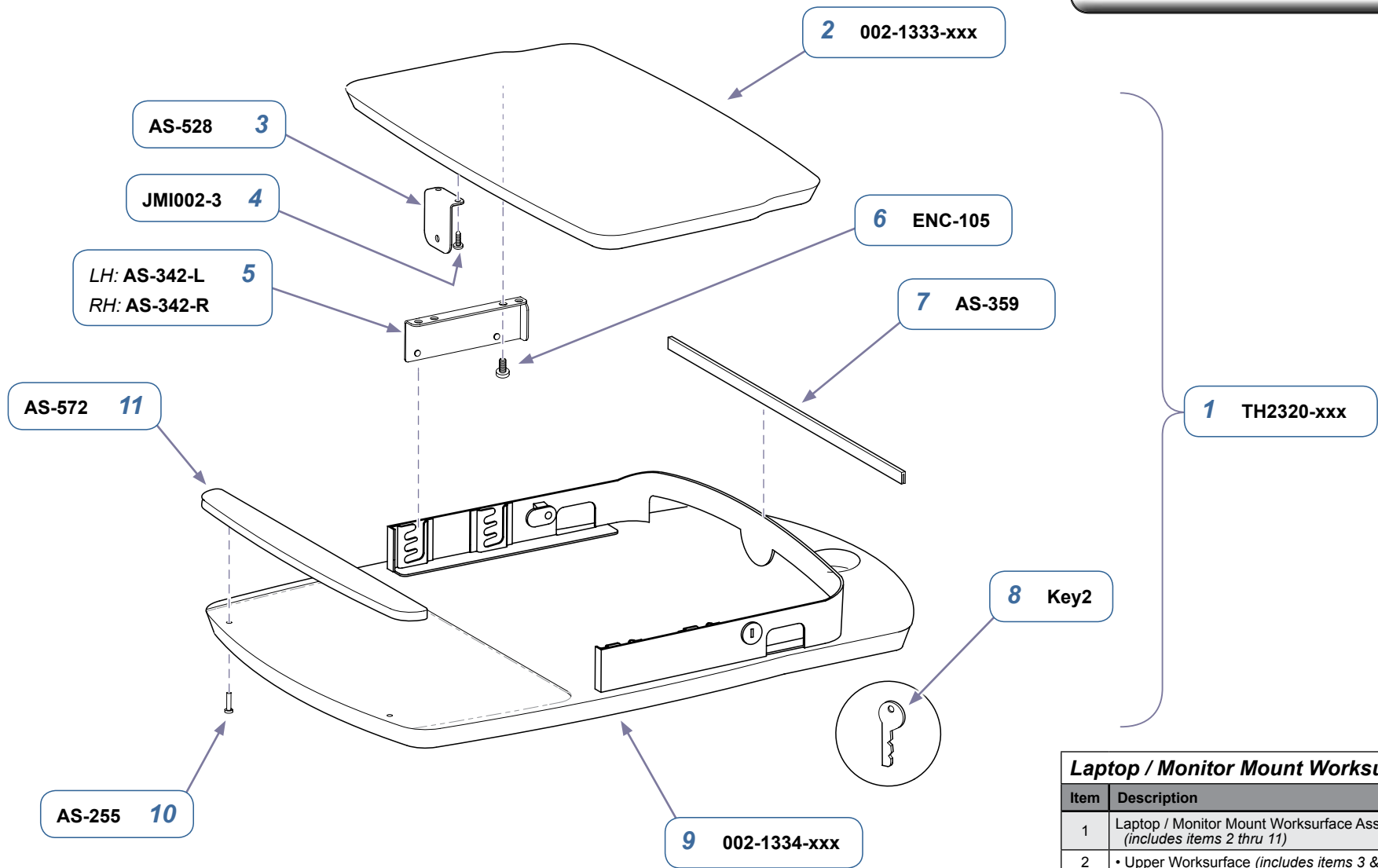


Laptop / Monitor Mount Worksurface		
Item	Description	Qty.
1	Laptop / Monitor Mount Worksurface Assembly (includes items 2 thru 11)	1
2	• Upper Worksurface (includes items 3 & 4)	1
3	• • Cord Retainer	1
4	• • Screw (#10 x 3/4")	2
5	• Upper Mounting Bracket (LH shown)	2
6	• Screw (1/4-20 x 5/8")	4
7	• Vinyl Edge Trim (per foot)	2 ft
8	• Key	1
9	• Lower Worksurface / Lock Assembly (includes items 10 & 11)	1
10	• • Screw (#8-32 x 1")	2
11	• • Wrist Rest	1

Models:	6203	
Serial Numbers:	<i>all</i>	

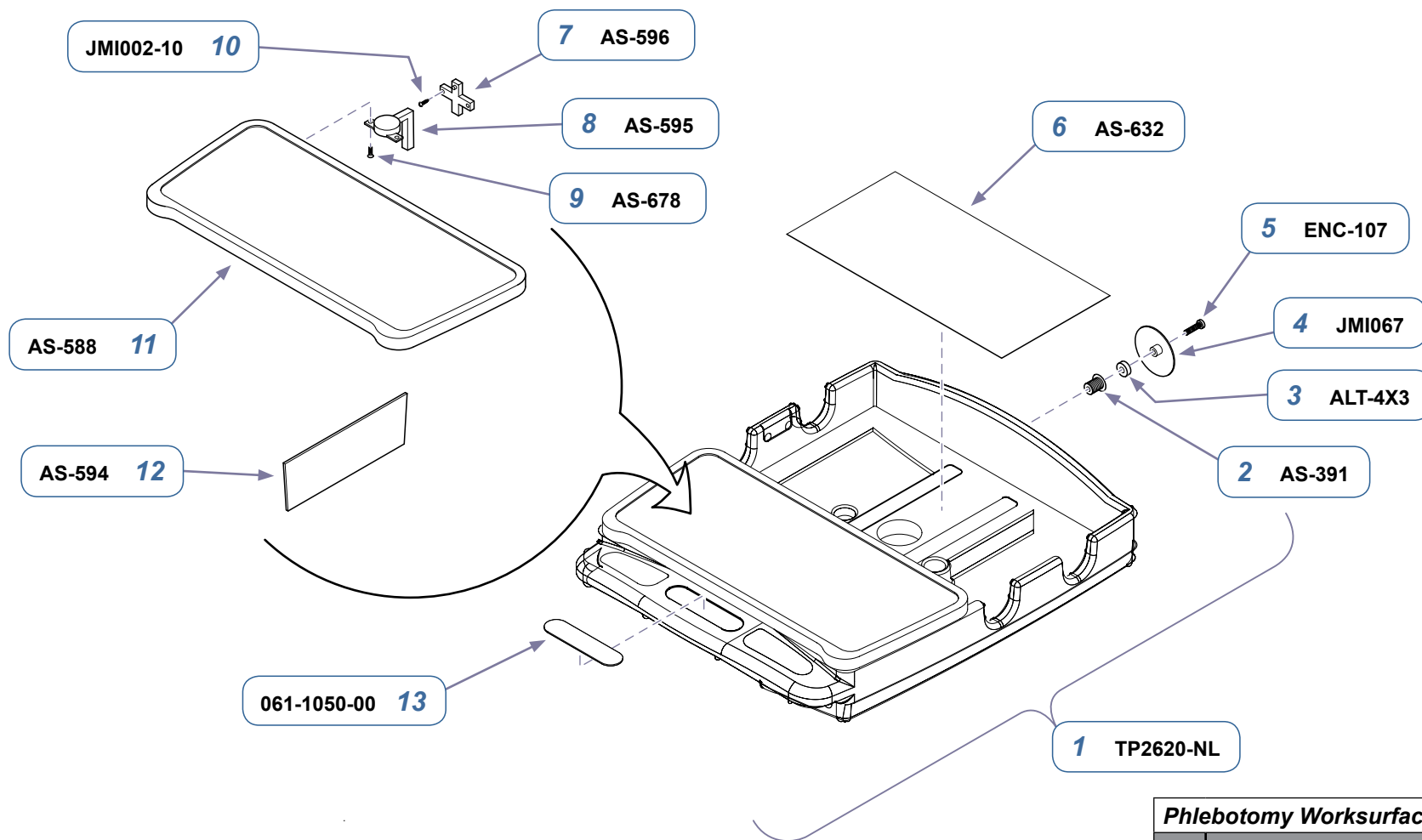



Phlebotomy Worksurface		
Item	Description	Qty.
1	Phlebotomy Worksurface Assembly (includes items 2 thru 13)	1
2	• Knife Thread Insert (1/4-20)	1
3	• Spacer	1
4	• Mouse Clamp	1
5	• Screw (1/4-20 x 1")	1
6	• Tub Insert	1
7	• Hinge Mounting Plate	2
8	• Top Hinge	2
9	• Screw (#8-32 x 1/2")	4
10	• Screw (#8 x 1/2")	4
11	• Flip-Up Worksurface	1
12	• Tub Divider	2
13	• Midmark Nameplate	1

Models:	6230	
Serial Numbers:	<i>all</i>	

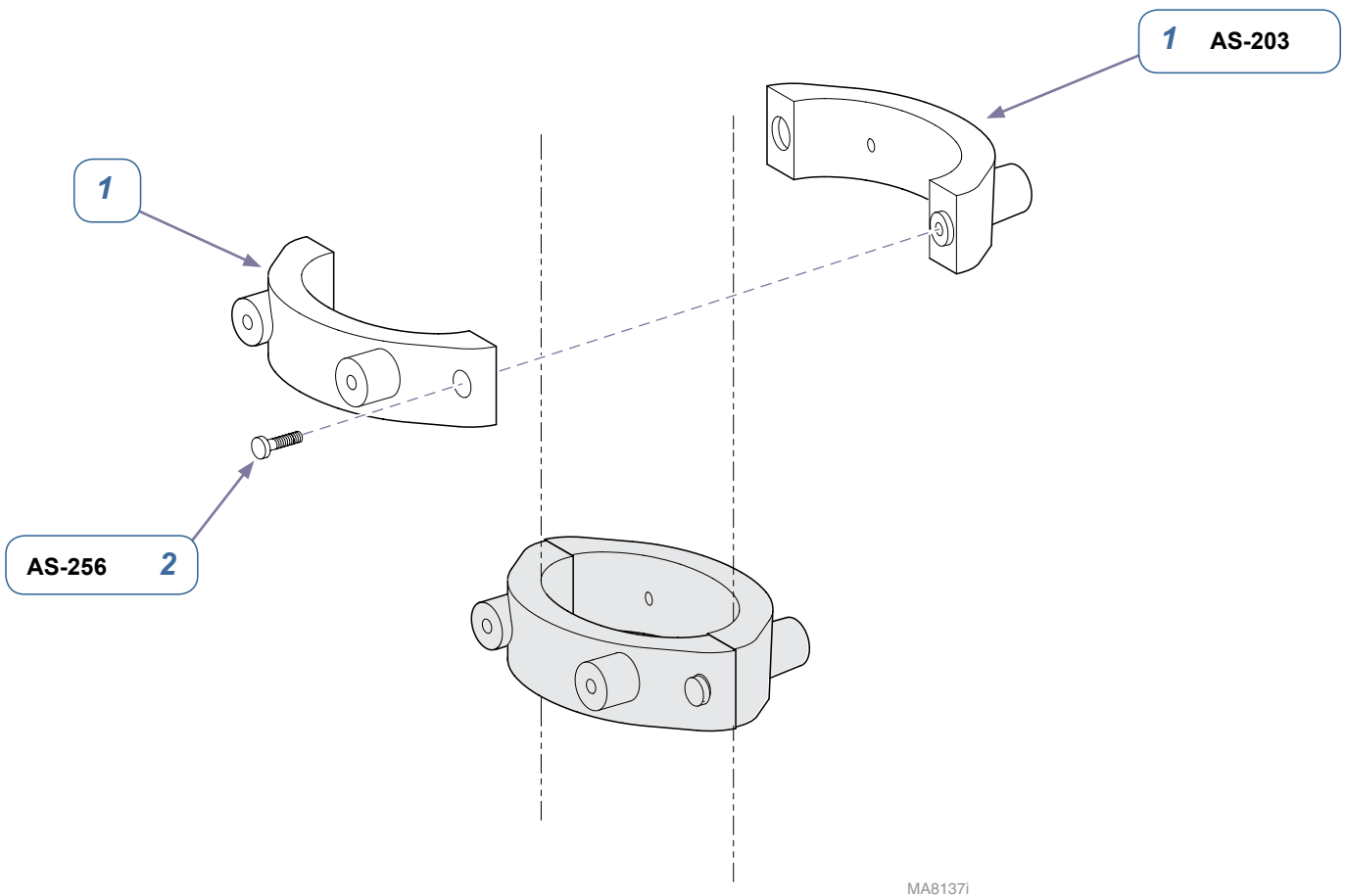

Models:	
Serial Numbers:	

Accessory Clamp Brackets		
Item	Description	Qty.
1	Accessory Clamp Bracket	4
2	Screw (#10-24 x 1")	4

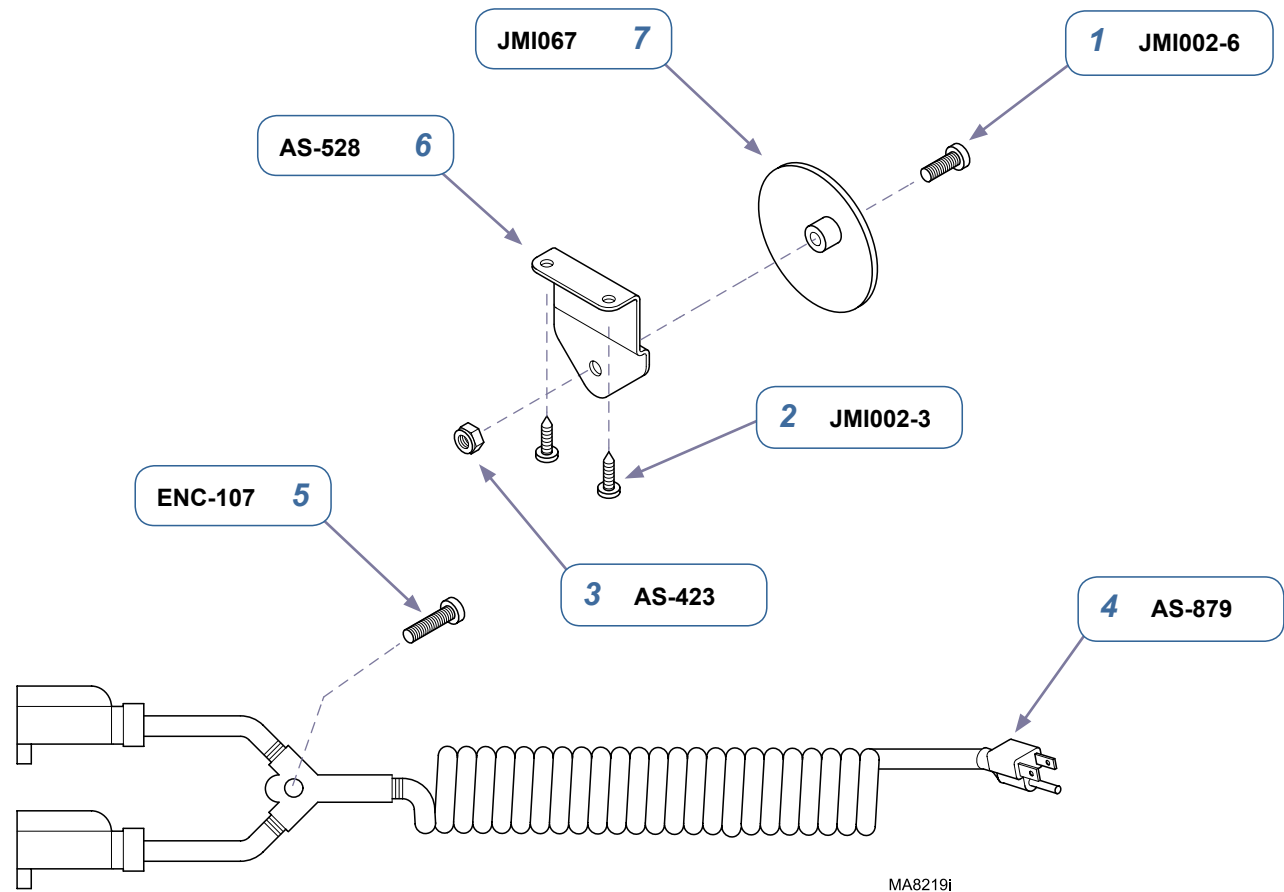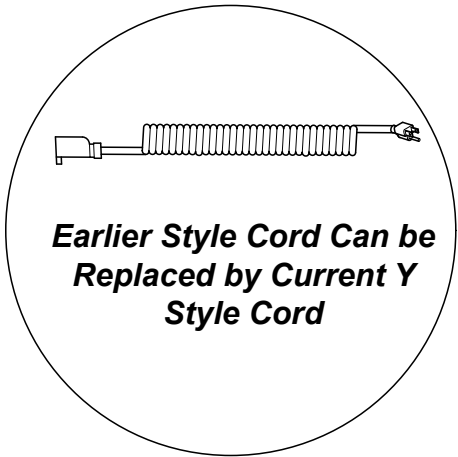


MA8219i

Coiled Power Cord Accessory - 9A435001 (includes items 1 thru 7)		
Item	Description	Qty.
1	Screw (1/4-20 x 3/4")	1
2	Screw (#10 x 3/4")	2
3	Hex Nut (1/4-20)	1
4	Coiled Power Cord	1
5	Screw (1/4-20 x 1")	1
6	Cord Retainer Bracket	1
7	Mouse Clamp	1

Models:

Serial Numbers:


Note: Retractable power cord is 6.5 ft. long.

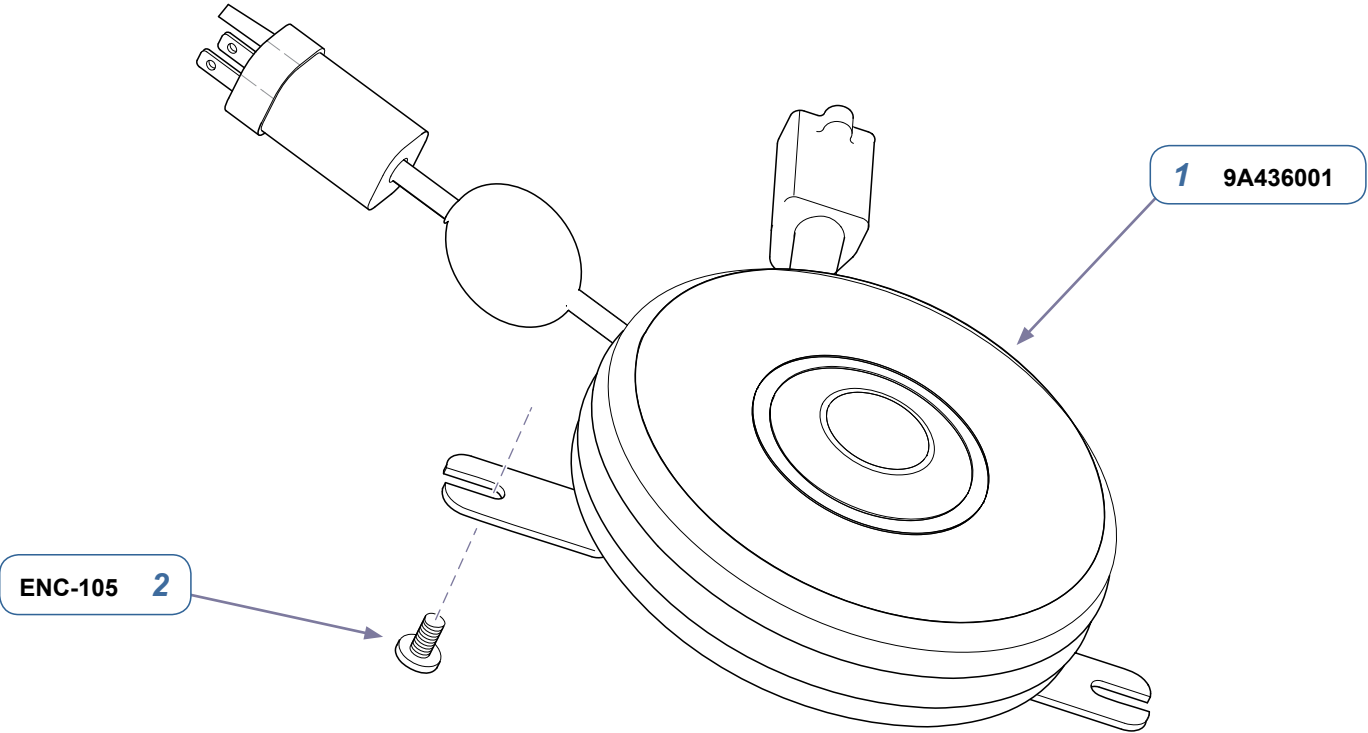

MA8145i

Models:	
Serial Numbers:	



Retractable Cable Reel - 9A436001		
Item	Description	Qty.
1	Retractable Cable Reel (includes item 2)	1
2	• Screw (1/4-20 x 5/8")	2

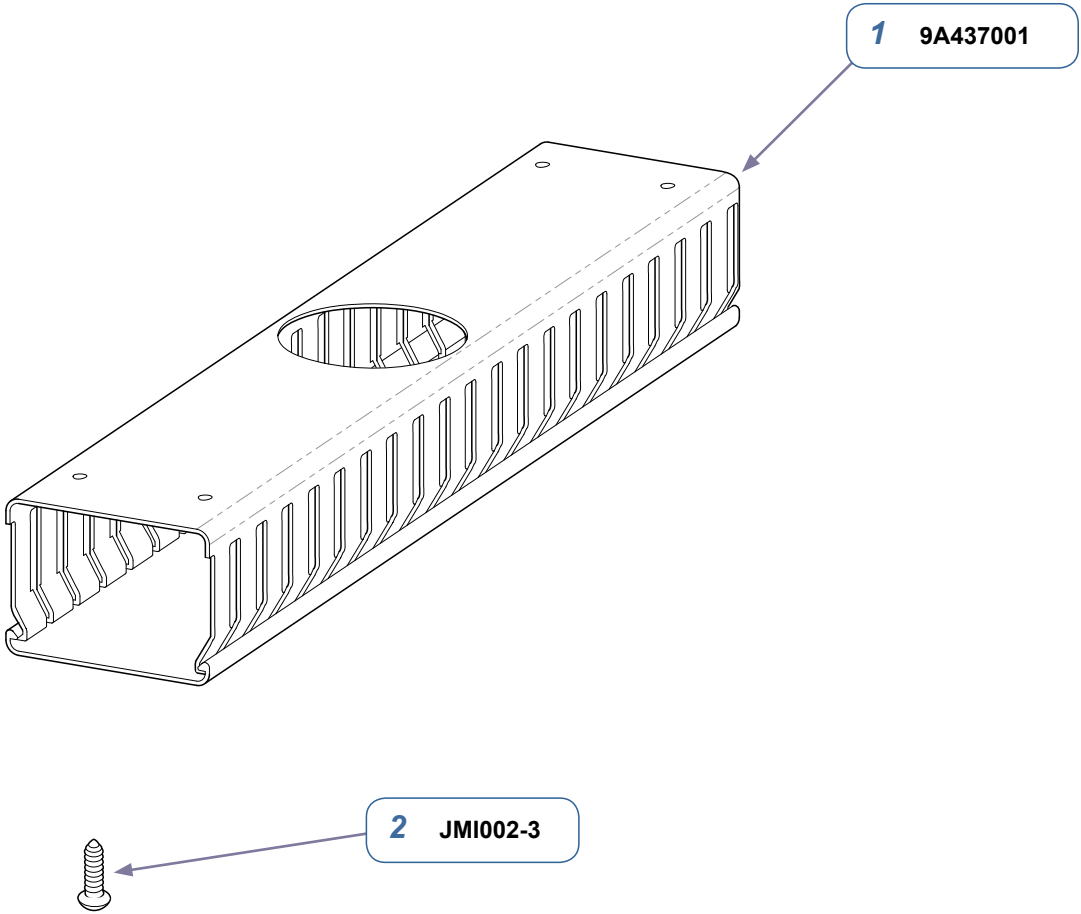

MA8147i

Models:	
Serial Numbers:	


Cable Management Holder - 9A437001		
Item	Description	Qty.
1	Cable Management Holder (includes item 2)	1
2	Screw (#10 x 3/4")	4

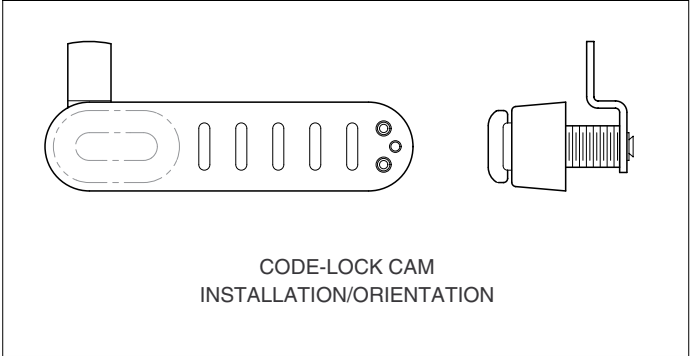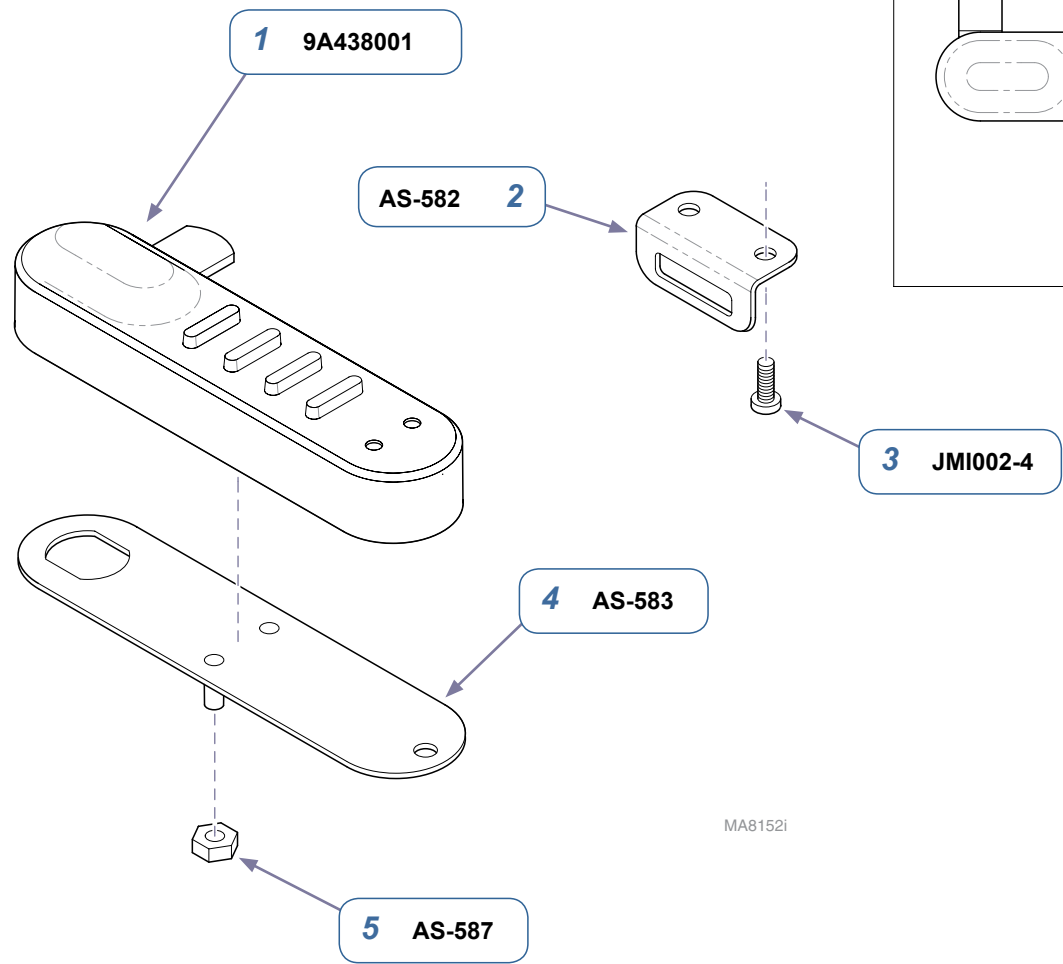

Models:

Serial Numbers:


Keyless Security System - 9A438001		
Item	Description	Qty.
1	Keyless Security System Accessory <i>(includes items 2 thru 5)</i>	
2	• Lock Bracket	1
3	• Screw (#8-32 x 1/2")	2
4	• Plate	1
5	• Nylock Hex Nut (#10-32)	2

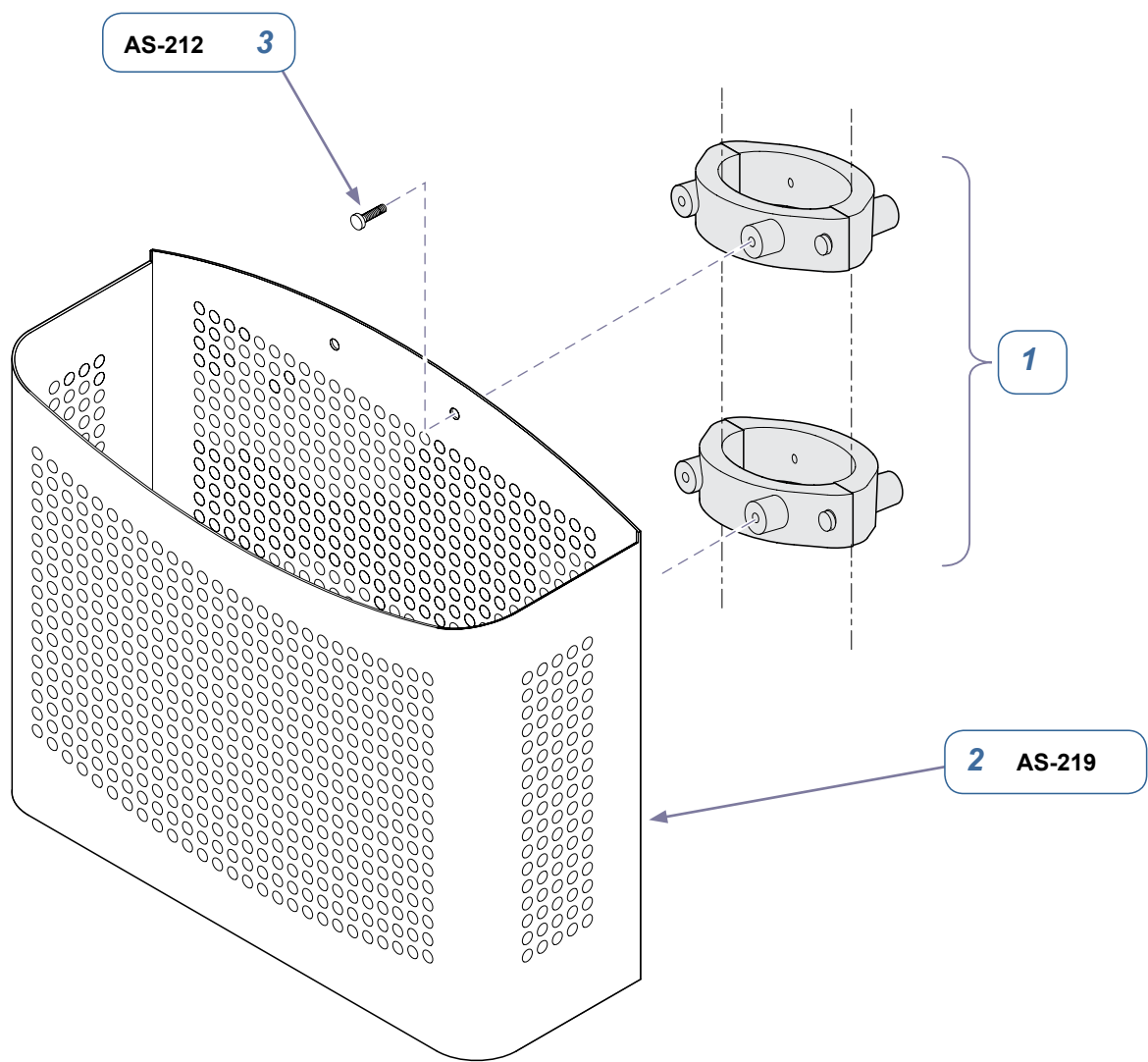

MA81411

Models:	
Serial Numbers:	

Utility Basket Accessory - 9A439001 (includes items 1 thru 3)		
Item	Description	Qty.
1	Refer to: " Accessory Mounting Brackets "	Ref
2	Utility Basket	1
3	Screw (#10-24 x 7/8")	4

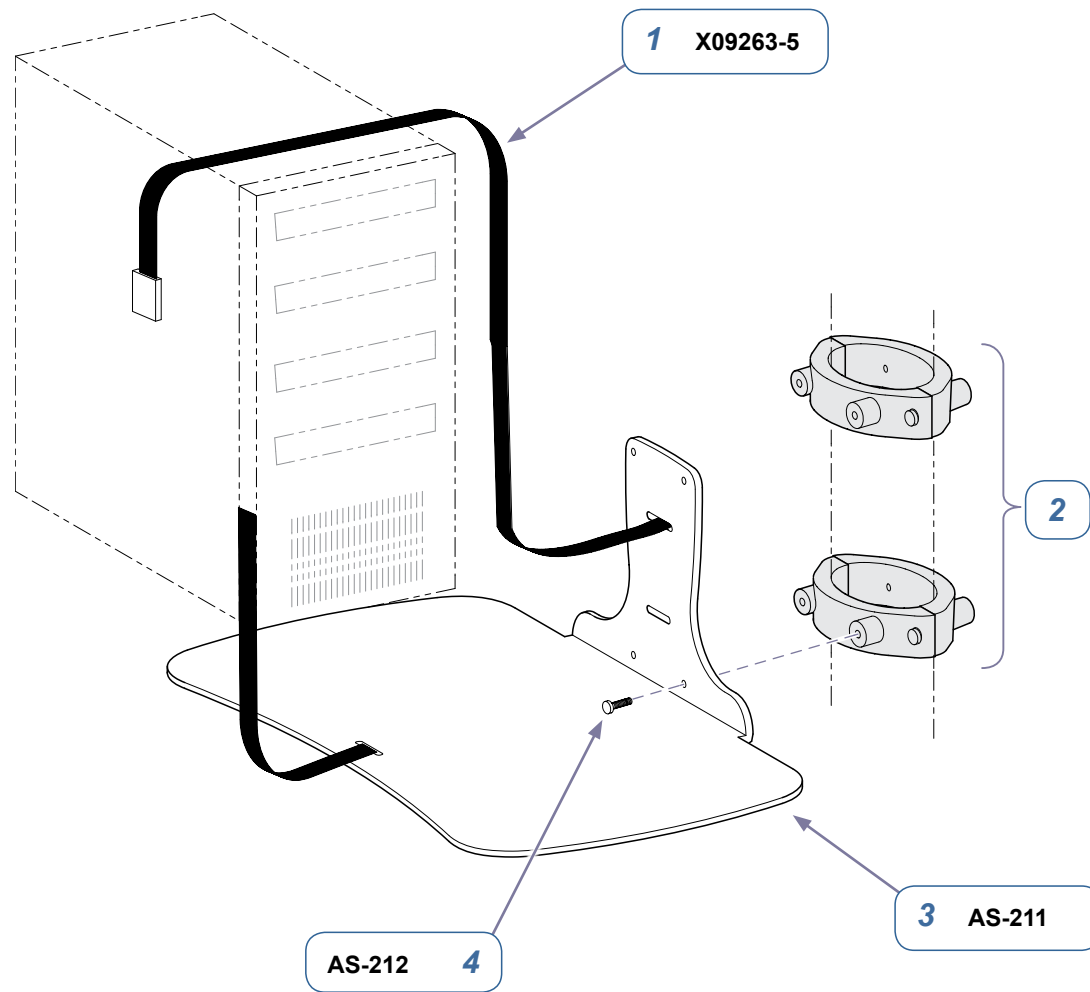


Shelf Dimensions: 20" w x 12" D

MA8148i

Models:	
Serial Numbers:	


Multi-Purpose Shelf - 9A440001 (includes items 1 thru 4)		
Item	Description	Qty.
1	Nylon Strap	1
2	Refer to: " Accessory Mounting Brackets "	Ref
3	Shelf	1
4	Hex Screw (#10-24 x 7/8")	4

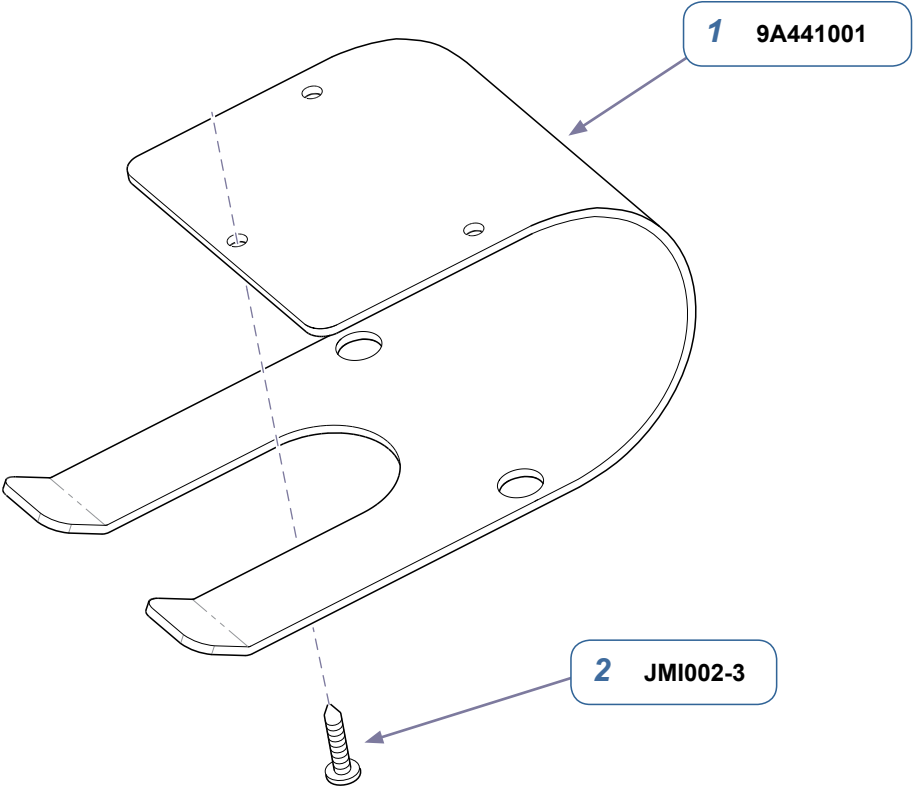

MA8162i


Models:	
Serial Numbers:	

Bar Code Holder - 9A441001		
Item	Description	Qty.
1	Bar Code Holder (includes item 2)	1
2	• Screw (#10 x 3/4")	3

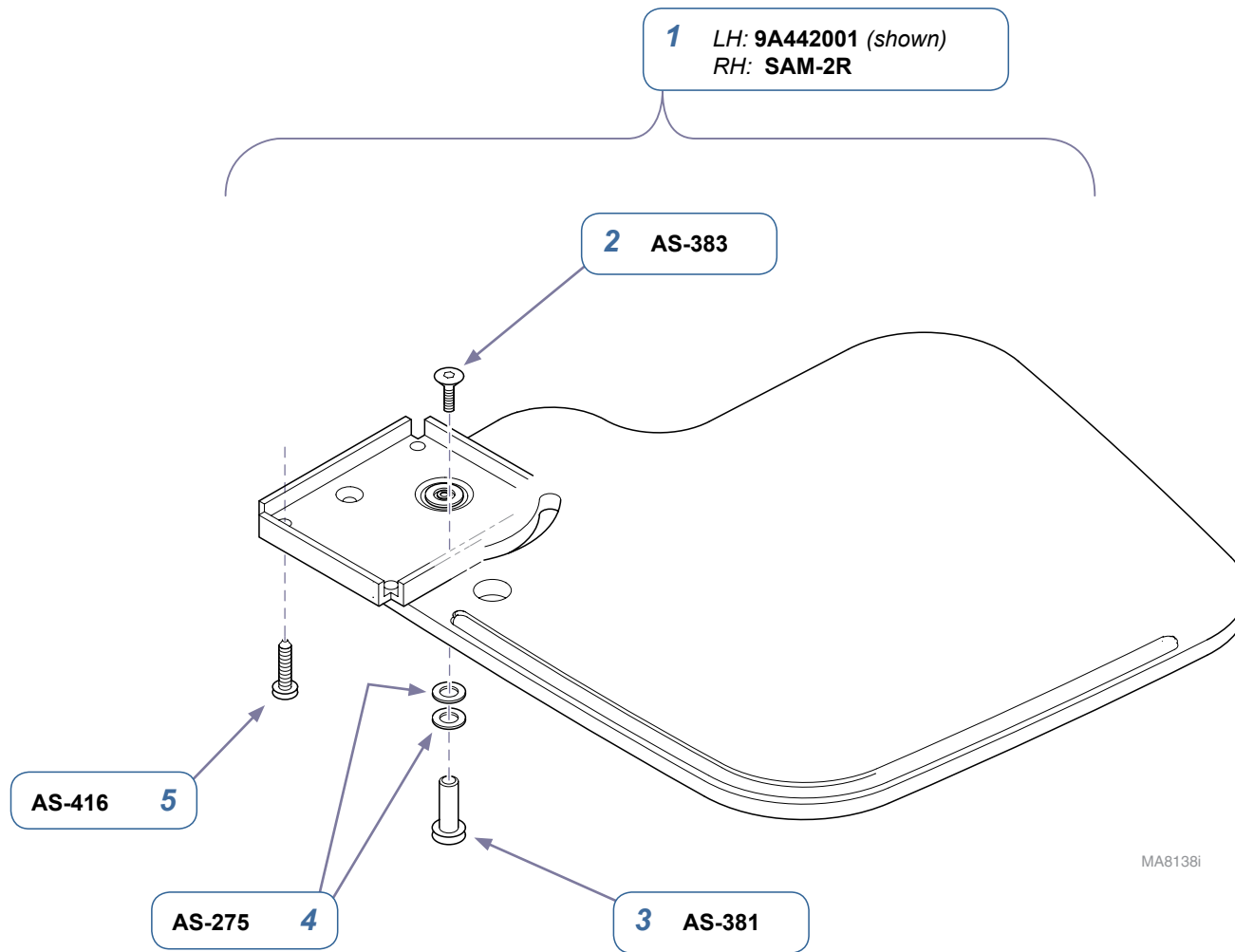


Swing Arm Mouse Pad		
Item	Description	Qty.
1	Swing Arm Mouse Pad Assy (incl. items 2 thru 5)	1
2	• Self-Locking Cap Screw (#10-24 x 1/2")	4
3	• Binding Barrel (#10-24 x 3/4")	1
4	• Nylon Pivot Washer	5
5	• Screw (#10 x 1")	1

Models:	all	
Serial Numbers:	prior to 4/1/2011	

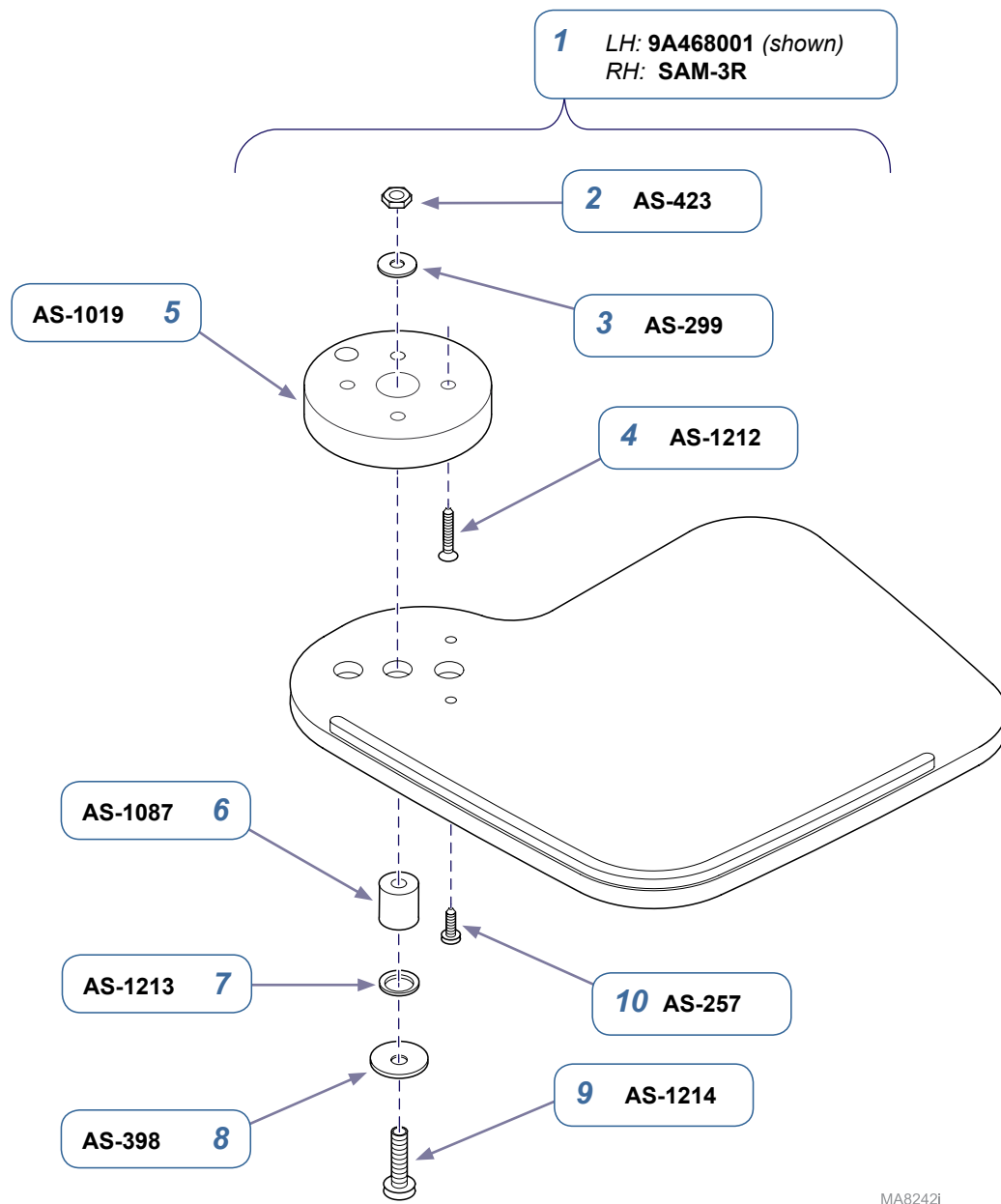


Swing Arm Mouse Pad		
Item	Description	Qty.
1	Swing Arm Mouse Pad Assy (incl. items 2 thru 10)	1
2	• Nylock Hex Nut (1/4-20)	1
3	• Washer (1/4")	1
4	• Screw (#10-24 x 1 1/4")	4
5	• Pivot Puck	1
6	• Pivot Bushing/Spacer	1
7	• Washer	1
8	• Washer	1
9	• Screw (1/4-20 x 1")	1
10	• Screw (#10-24 x 1/2")	2

Models:	all	
Serial Numbers:	after 4/1/2011	

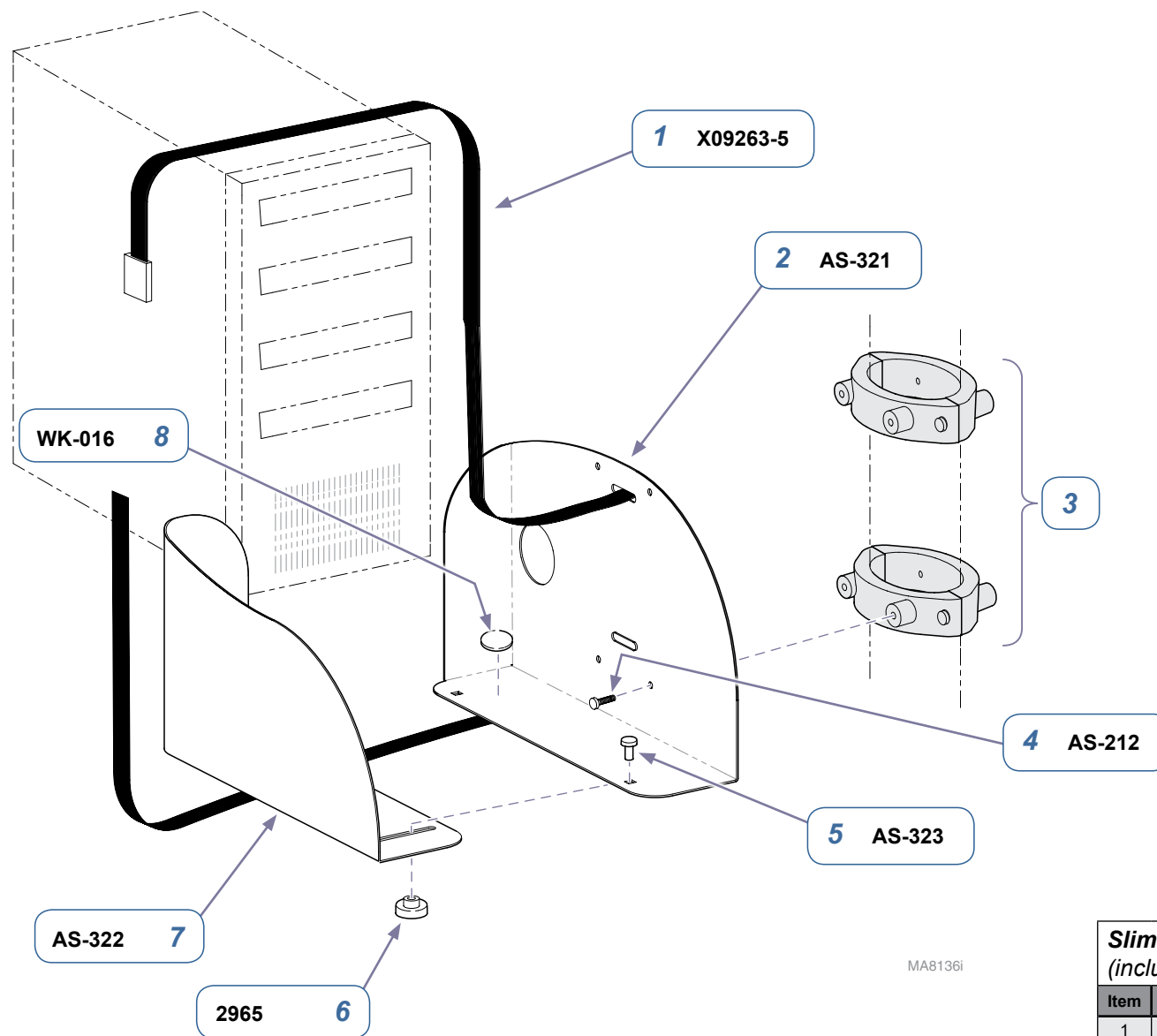


Slim Line CPU Holder - 9A443001 (includes items 1 thru 8)		
Item	Description	Qty.
1	Nylon Strap	1
2	CPU Holder - Main Body	1
3	Refer to: " Accessory Mounting Brackets "	Ref
4	Screw (#10-24 x 7/8")	4
5	Carriage Bolt (1/4-20 x 1/2")	2
6	Knob (1/4-20 thread)	2
7	CPU Holder - Clamp	1
8	Non-Skid Pad	2

Models:	all	
Serial Numbers:	prior to 12/12/2011	

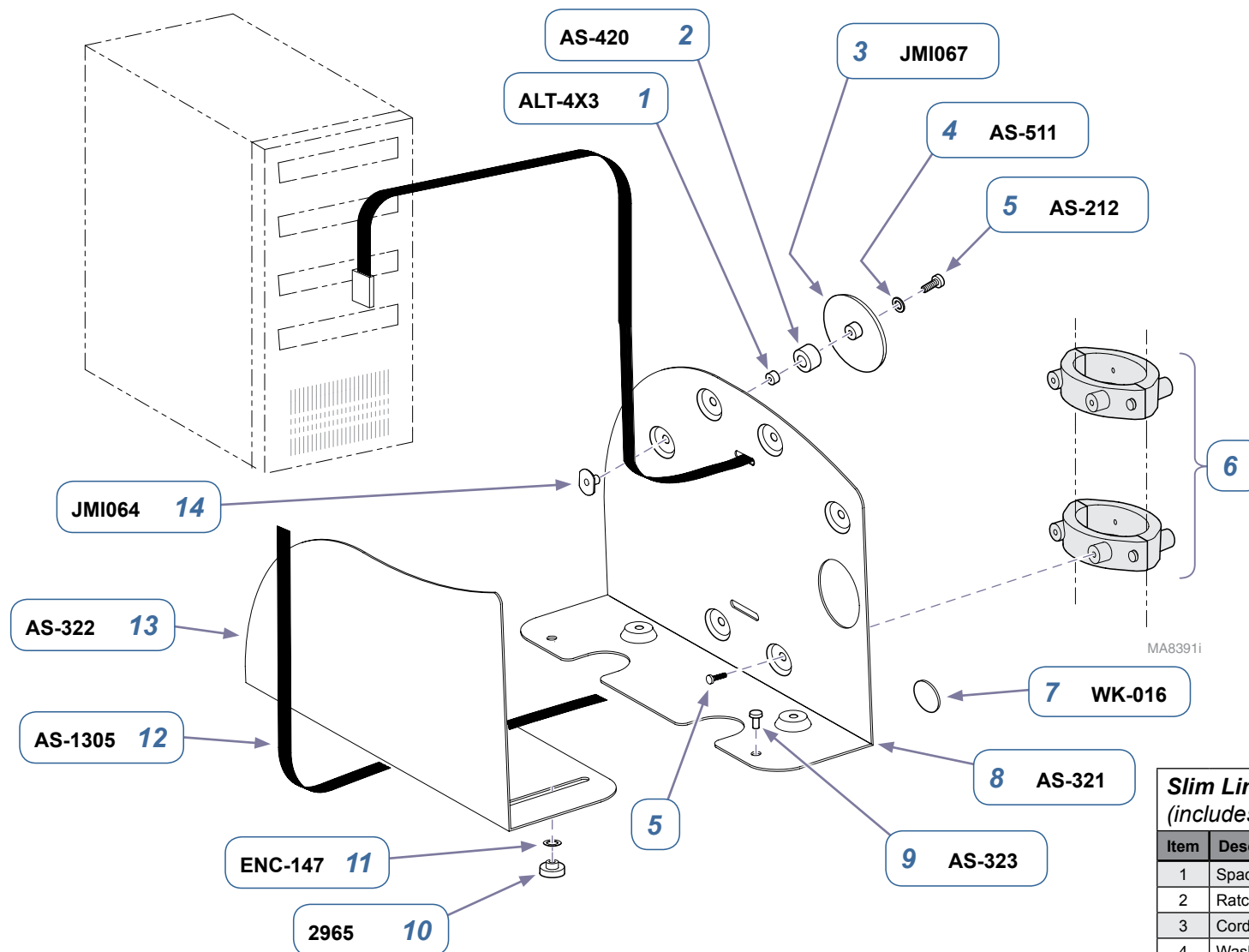

Slim Line CPU Holder - 9A443001 (includes items 1 thru 8)		
Item	Description	Qty.
1	Spacer	2
2	Ratchet Wheel Spacer	2
3	Cord Puck	2
4	Washer (#10)	2
5	Screw (#10-24 x 7/8")	6
6	Refer to: " Accessory Mounting Brackets "	Ref
7	Non-Skid Pad	8
8	CPU - Main Body	1
9	Carriage Bolt (1/4-20 x 1/2")	2
10	Knob	2
11	Washer (5/16")	2
12	Nylon Strap	1
13	CPU - Clamp	1
14	Weldnut	2


Models:	all	
Serial Numbers:	after 12/12/2011	

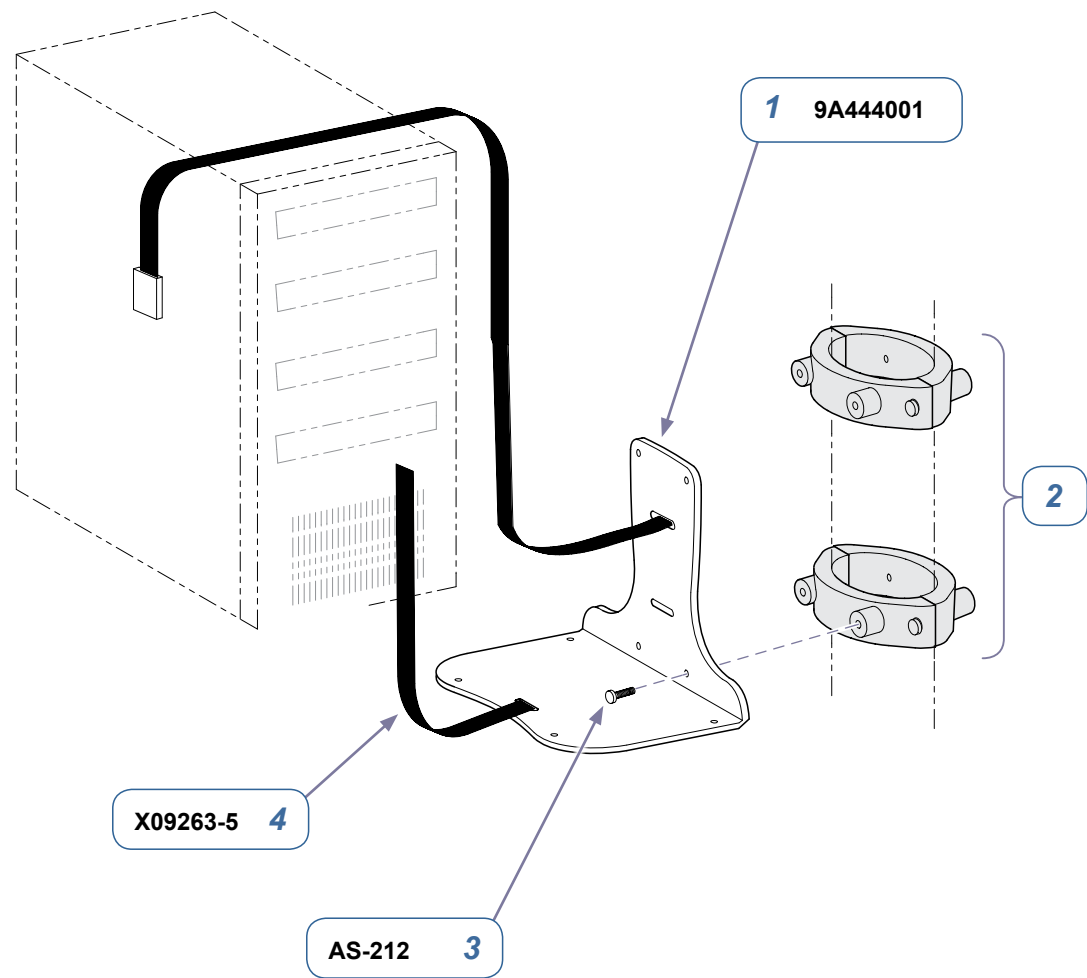

MA8149i

Models:	
Serial Numbers:	


CPU Holder Accessory - 9A444001		
Item	Description	Qty.
1	CPU Holder Assembly (includes items 2 thru 4)	1
2	• Refer to: " Accessory Mounting Brackets "	Ref
3	• Screw (#10-24 x 7/8")	4
4	• Nylon Strap	1

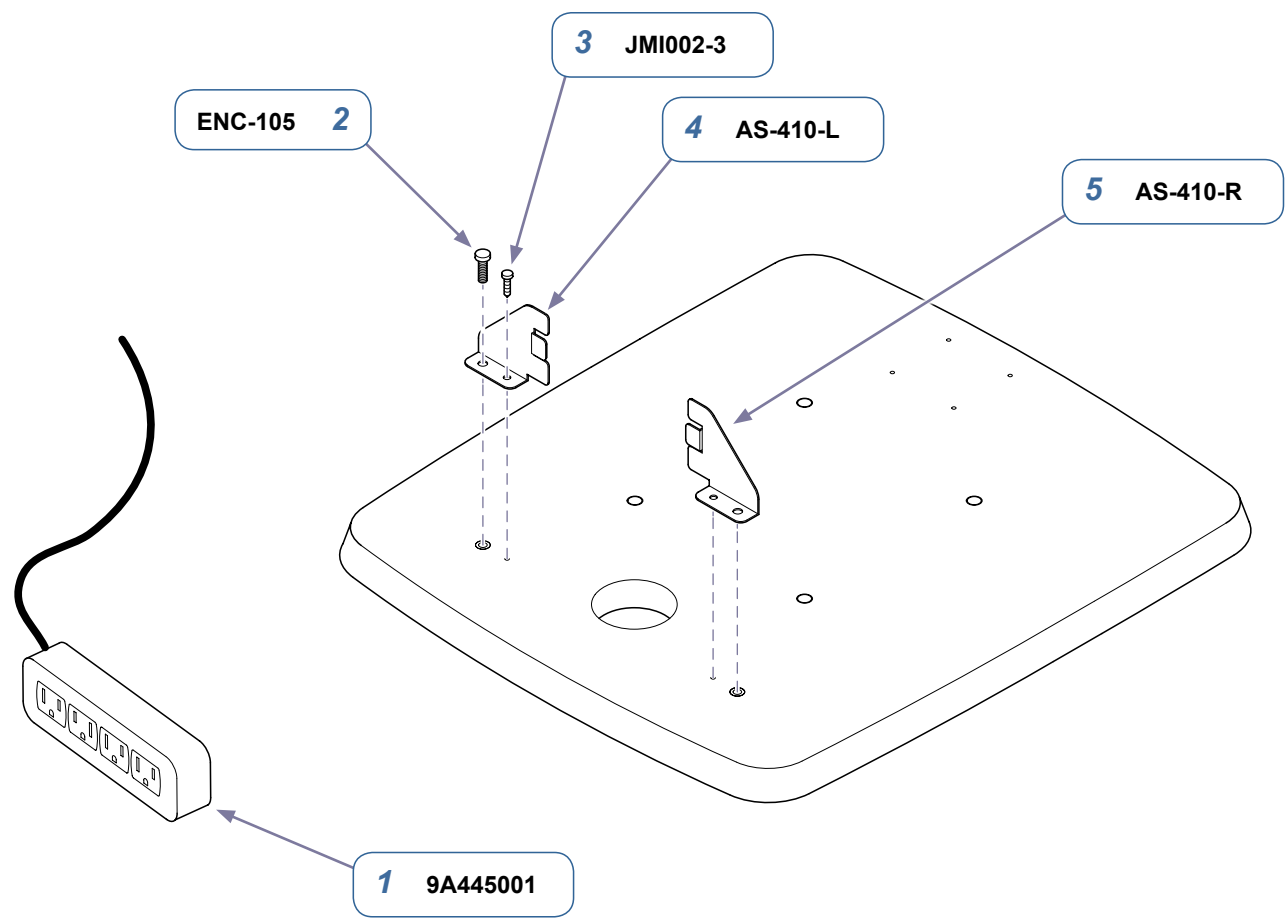

MA8151i

Models:	
Serial Numbers:	


Power Distribution Strip Accessory - 9A445001 (includes items 1 thru 5)		
Item	Description	Qty.
1	Power Distribution Strip Accy (incl. items 2 thru 5)	1
2	• Pan Head Screw (1/4-20 x 5/8")	2
3	• Screw (#10 x 3/4")	2
4	• Mounting Bracket - Left	1
5	• Mounting Bracket - Right	1

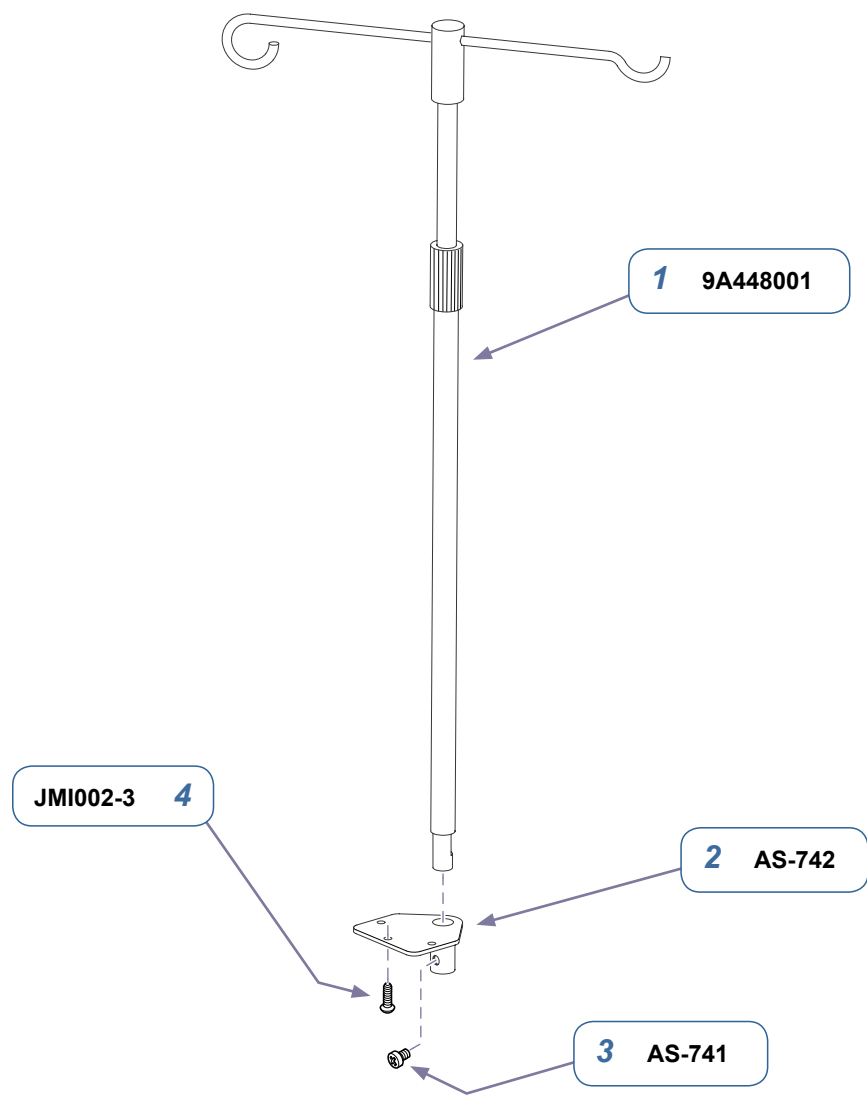

MA8150i

Models:	
Serial Numbers:	


IV Pole / ECG Holder - 9A448001		
Item	Description	Qty.
1	IV Pole / ECG Holder (includes items 2 thru 4)	1
2	• Bracket	1
3	• Screw (1/4-20 x 3/8")	1
4	• Screw (#10 x 3/4")	3

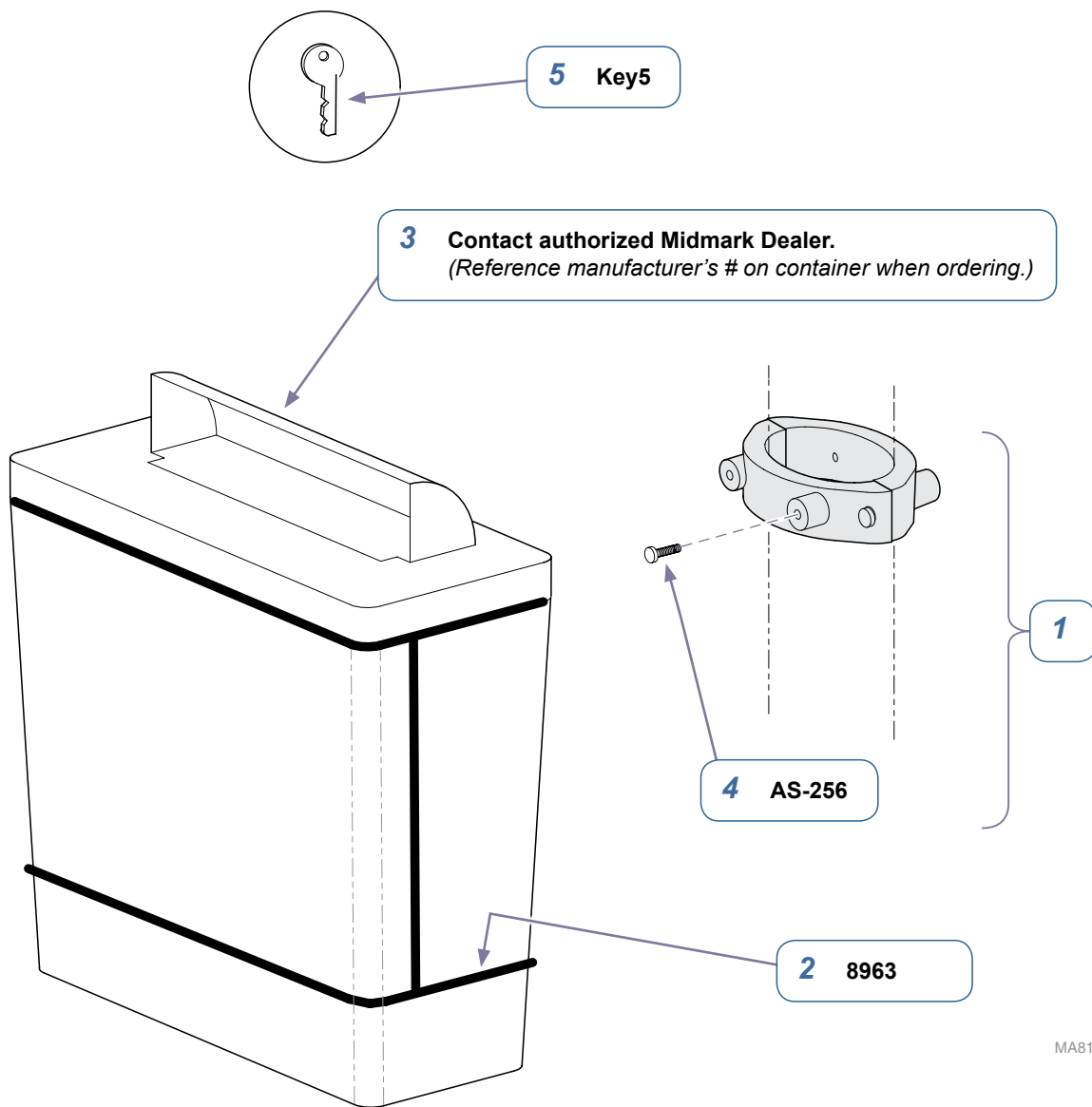

MA8142i

Models:	
Serial Numbers:	

Sharps Container Accessory - 9A450001 (includes items 1 thru 5)		
Item	Description	Qty.
1	Refer to: " Accessory Mounting Brackets "	Ref
2	Wire Hanger	1
3	Sharps Container	1
4	Screw (#10-24 x 7/8")	2
5	Key	2

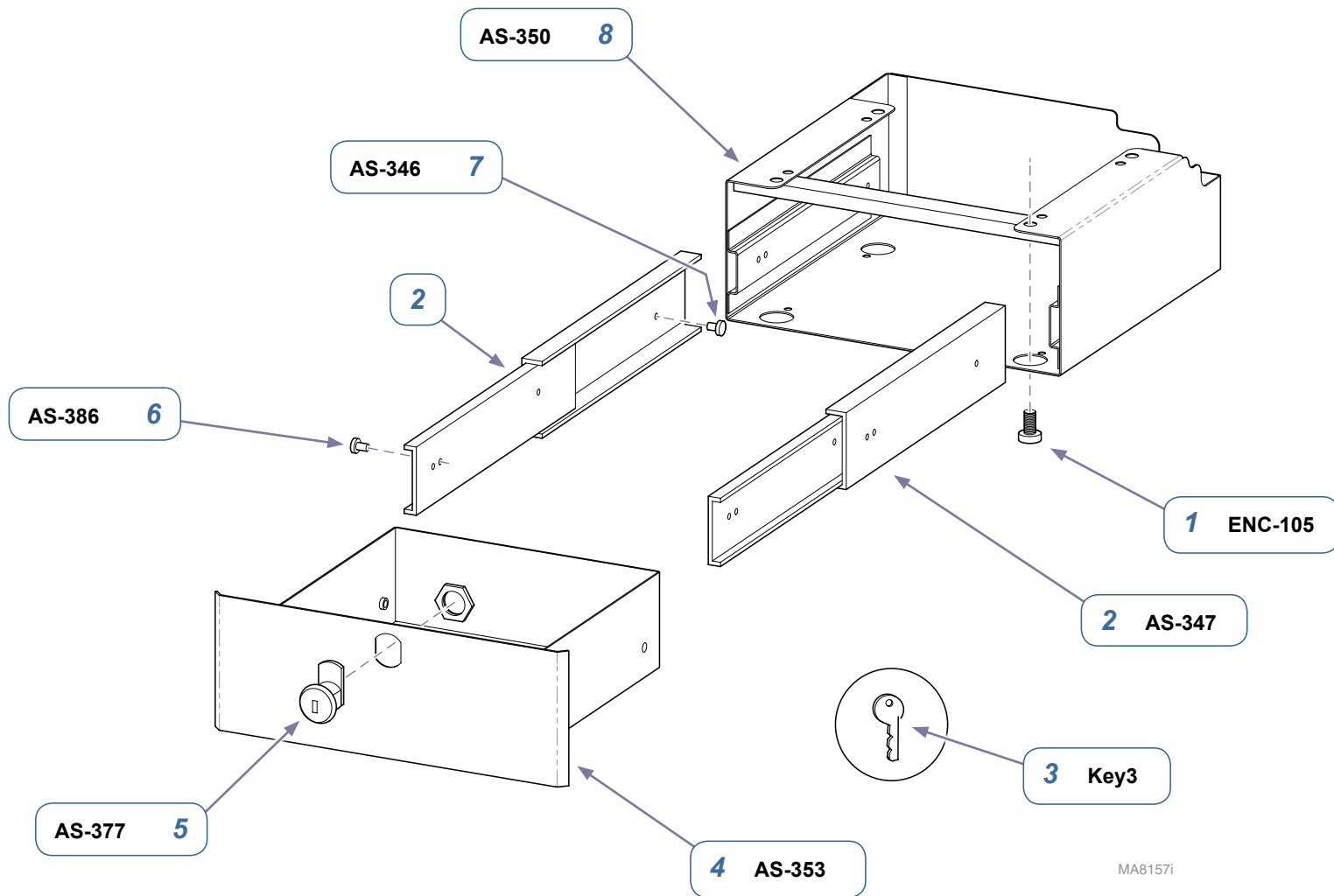


MA81571

Locking Drawer Accessory - 9A451001
(includes items 1 thru 8)

Item	Description	Qty.
1	Screw (1/4-20 x 5/8")	4
2	Drawer Slide (pair)	1
3	Key	1
4	Drawer Weldment	1
5	Lock Kit	1
6	Screw (#8-32 x 1/4")	4
7	Screw (#8 x 1/4")	4
8	Drawer Housing	1

Models:


Serial Numbers:

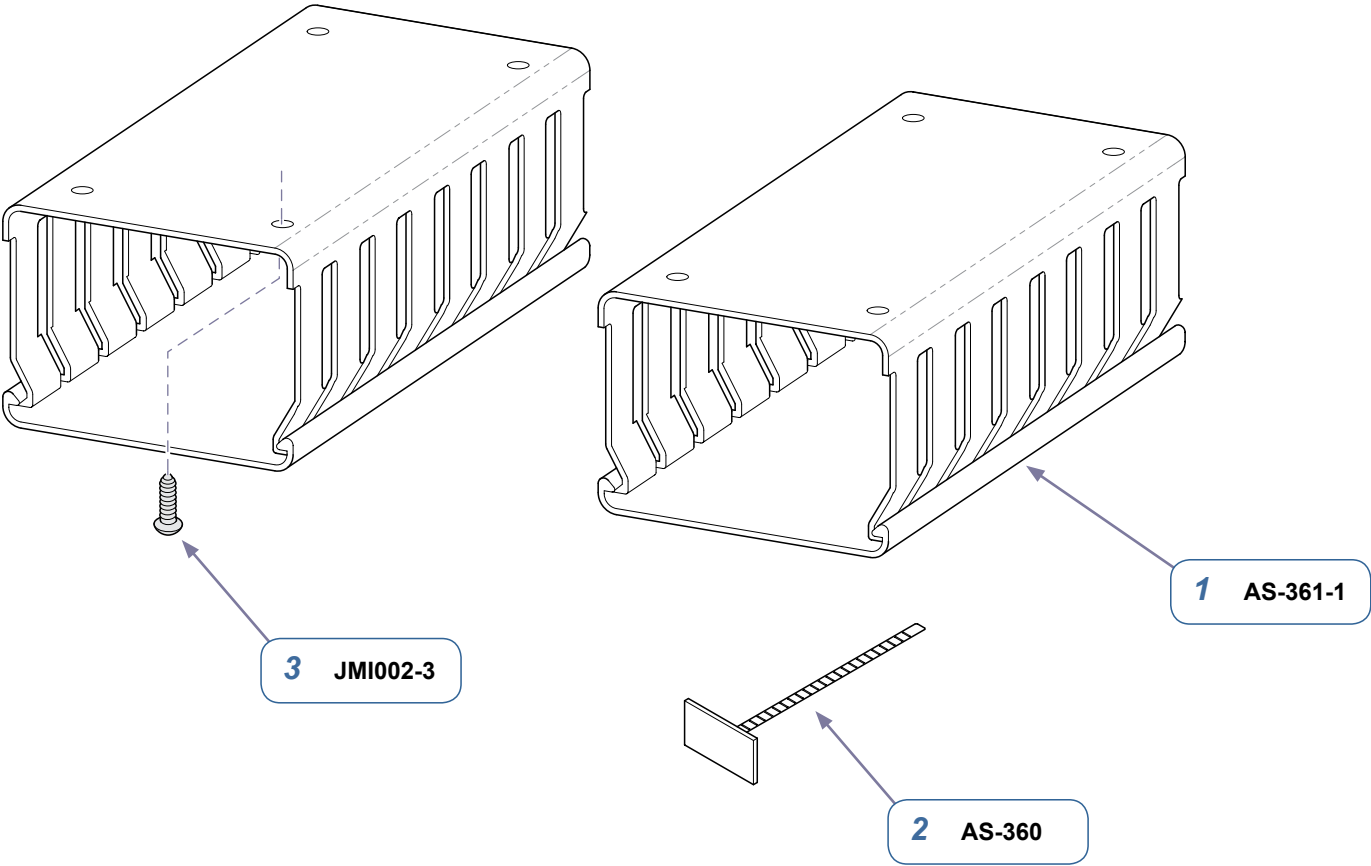

MA8166i

Models:	
Serial Numbers:	


Cable Management Accessory - 9A452001 (includes items 1 thru 3)		
Item	Description	Qty.
1	Wire Duct	2
2	Wire Tie (Adhesive Backed)	2
3	Screw (#10 x 3/4")	8

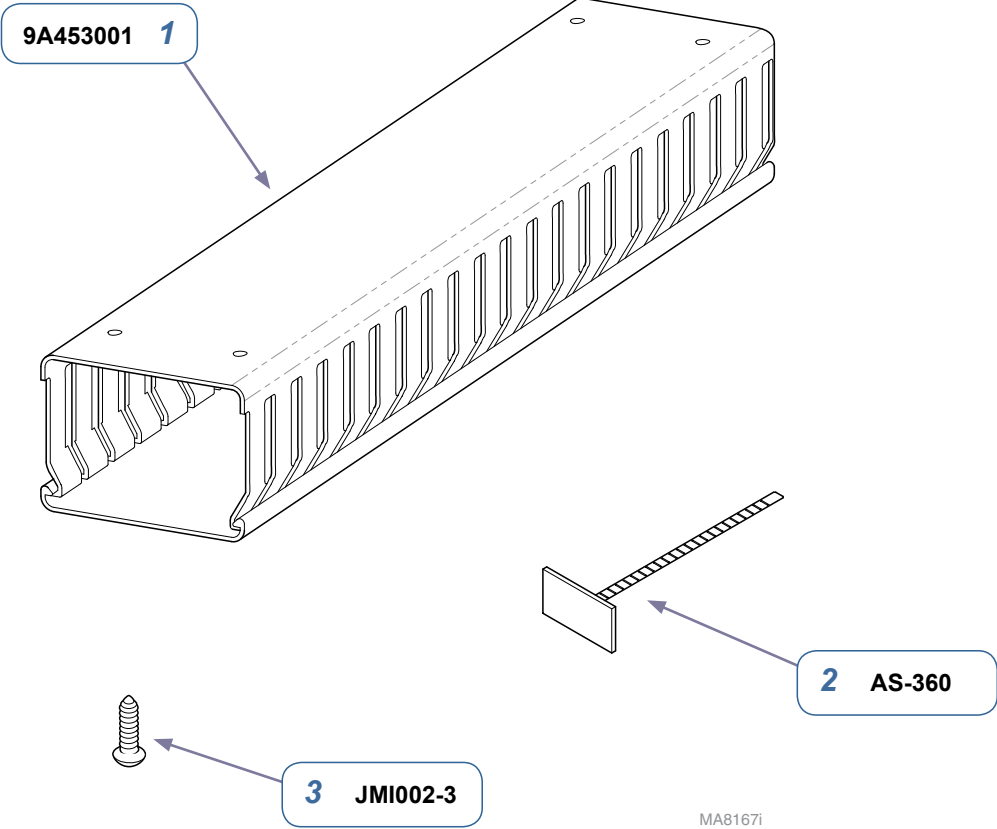


Models:	
Serial Numbers:	

Cable Management Accessory - 9A453001 (includes items 1 thru 3)		
Item	Description	Qty.
1	Wire Duct (includes items 2 & 3)	1
2	Wire Tie (Adhesive Backed)	2
3	Screw (#10 x 3/4")	8

Built before: 5/1/2011

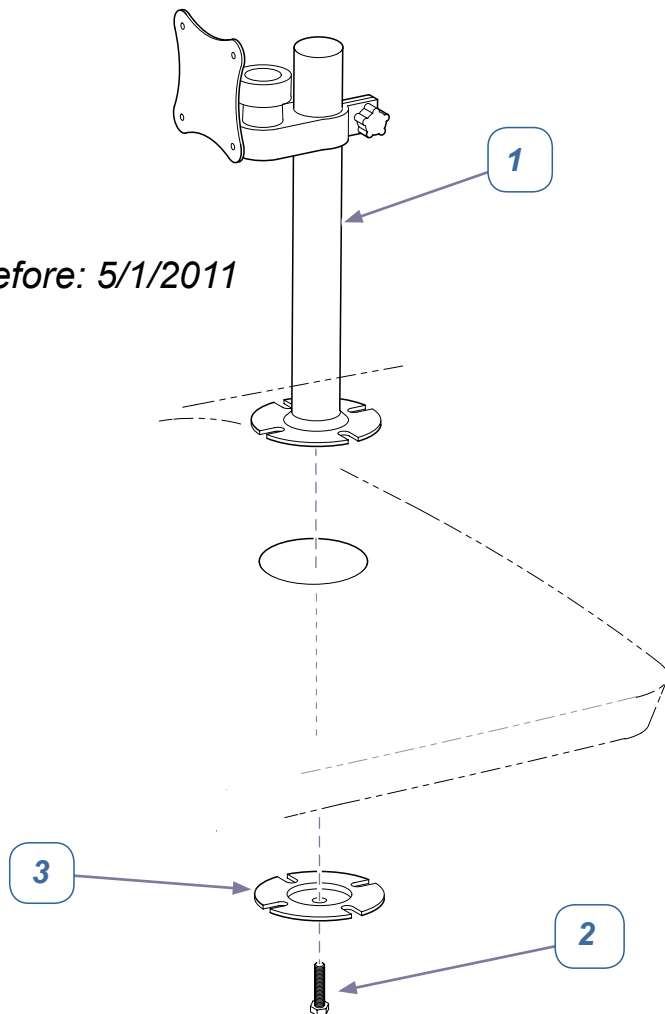


Built after: 5/1/2011

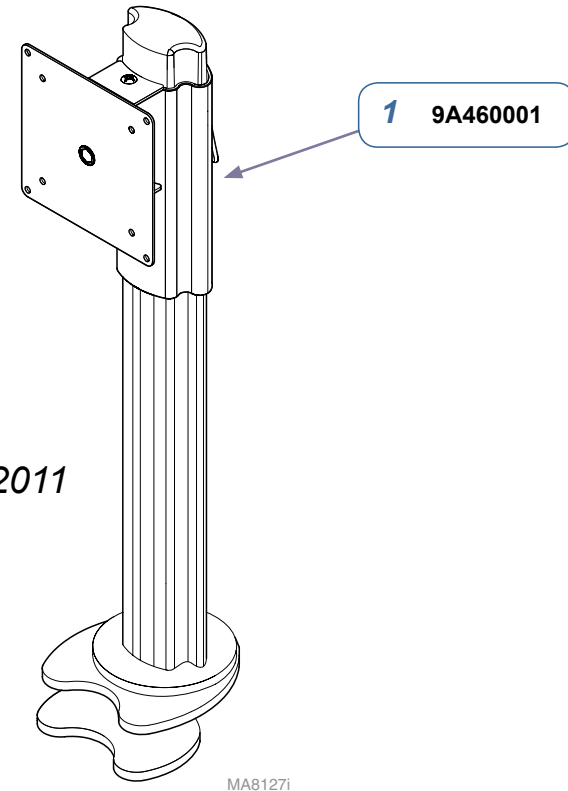

Models:	
Serial Numbers:	


Flat Panel Monitor Mount - 9A460001		
Item	Description	Qty.
1	Flat Panel Monitor Mount	
2	• Hex Bolt (3/8" x 3")	1
3	• Retainer Plate	1
4	• Wire Tie (not shown)	2

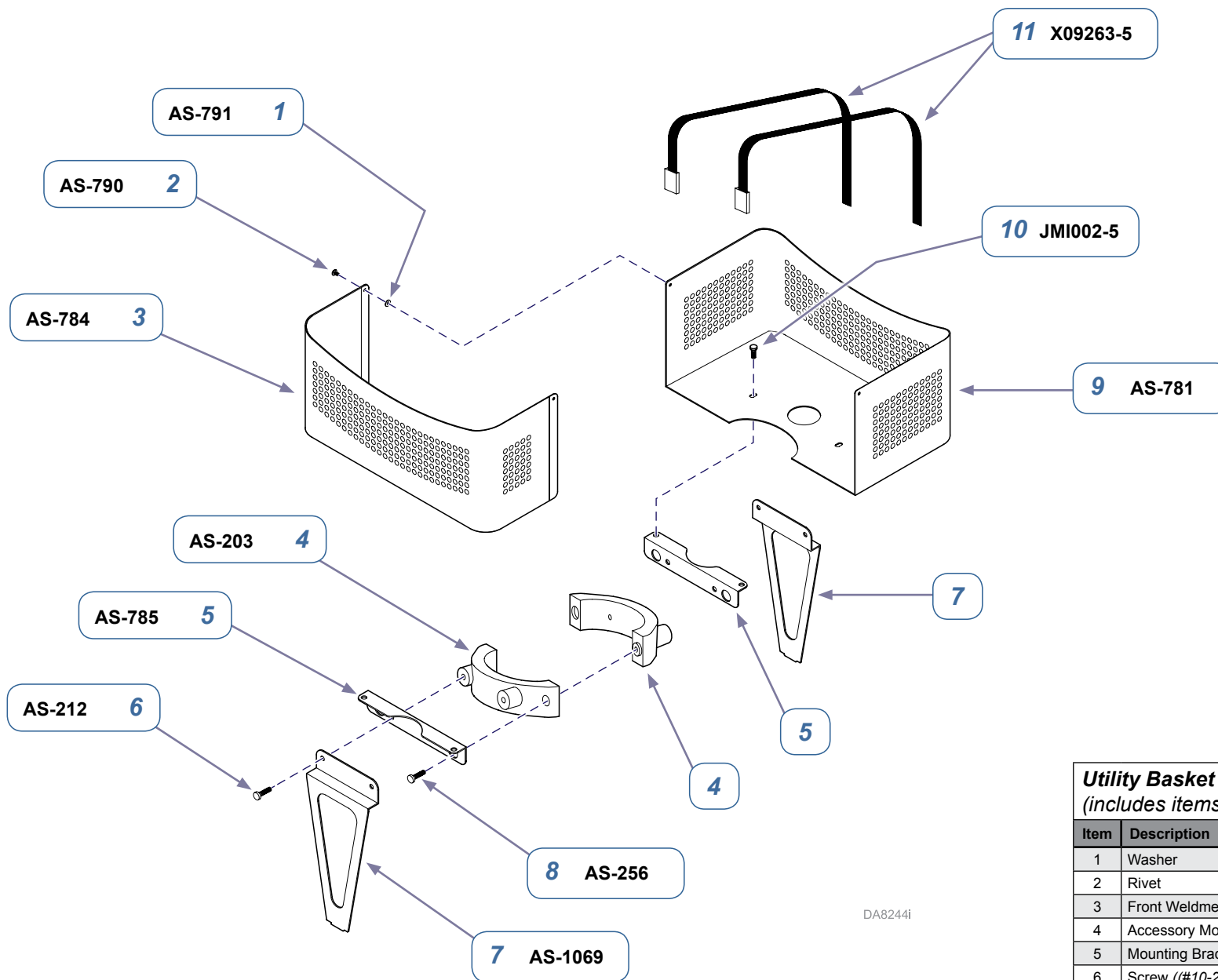

DA8244i

Utility Basket Accessory - 9A466001 (includes items 1 thru 11)		
Item	Description	Qty.
1	Washer	2
2	Rivet	2
3	Front Weldment	1
4	Accessory Mounting Bracket	2
5	Mounting Bracket	2
6	Screw (#10-24 x 7/8")	4
7	Spacer Bracket	2
8	Screw (#10-24 x 1")	2
9	Rear Weldment	1
10	Screw (#10-24 x 1/2")	4
11	Nylon Strap (24")	2

Models:	
Serial Numbers:	

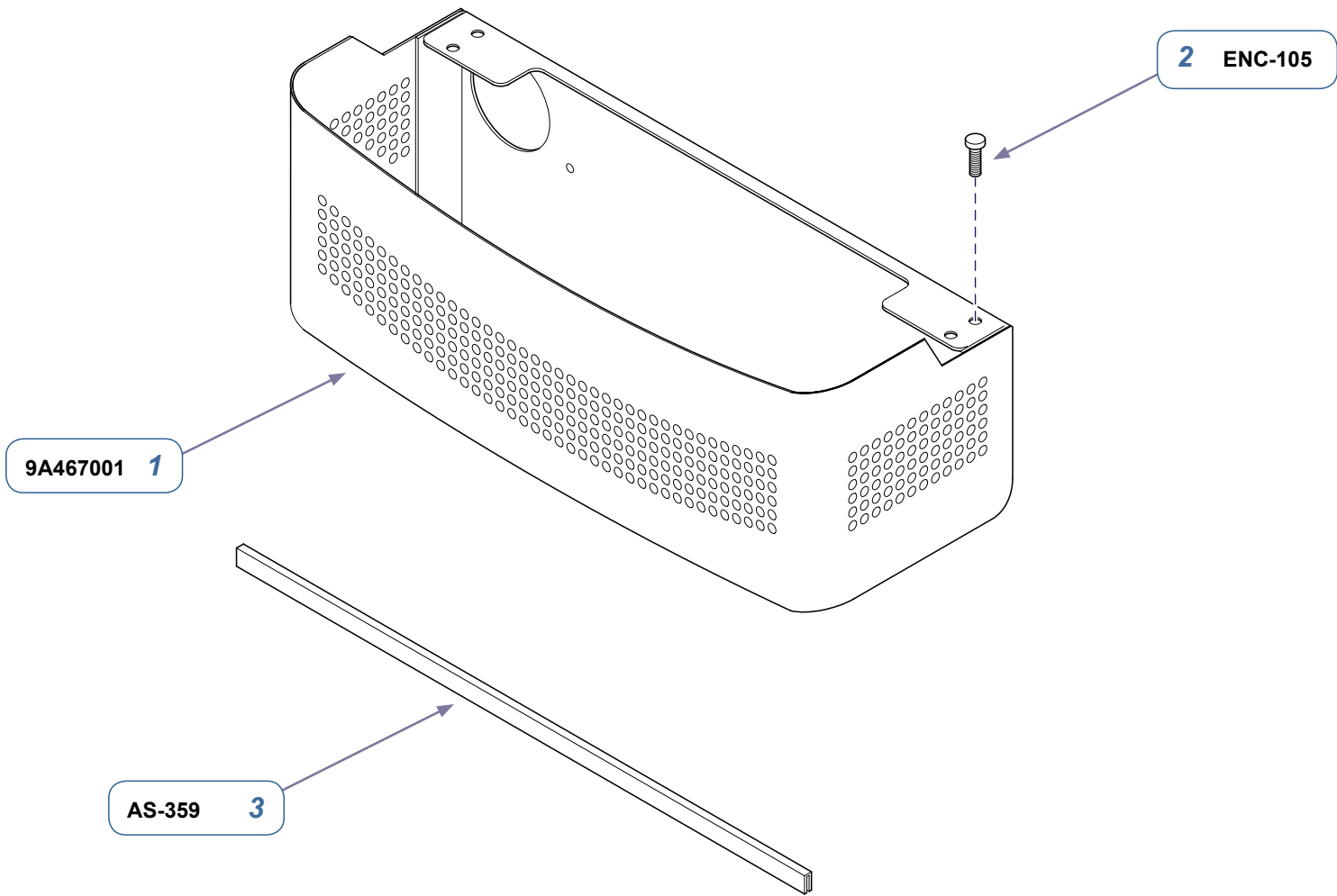


MA8243i


Models:	
Serial Numbers:	

Rear Basket Accessory - 9A467001		
Item	Description	Qty.
1	Rear Basket Assy. (includes items 2 & 3)	1
2	• Screw (#1/4-20 x 5/8")	4
3	• Vinyl Edge Trim (sold by the foot)	2


Drawer / Bin Options

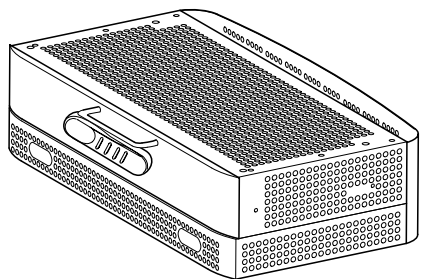

Med Drawer Storage
w/ Lock E-37

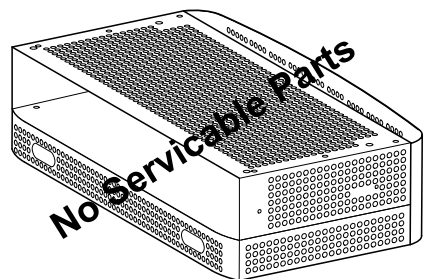


Med Drawer Chart Pocket

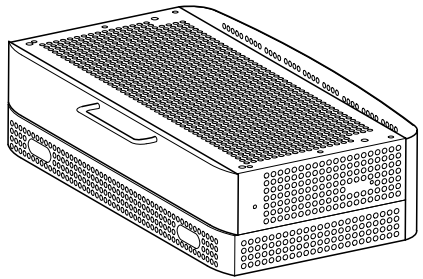

Med Drawer Storage
w/o Lock E-38


MA8248i

Models:	
Serial Numbers:	

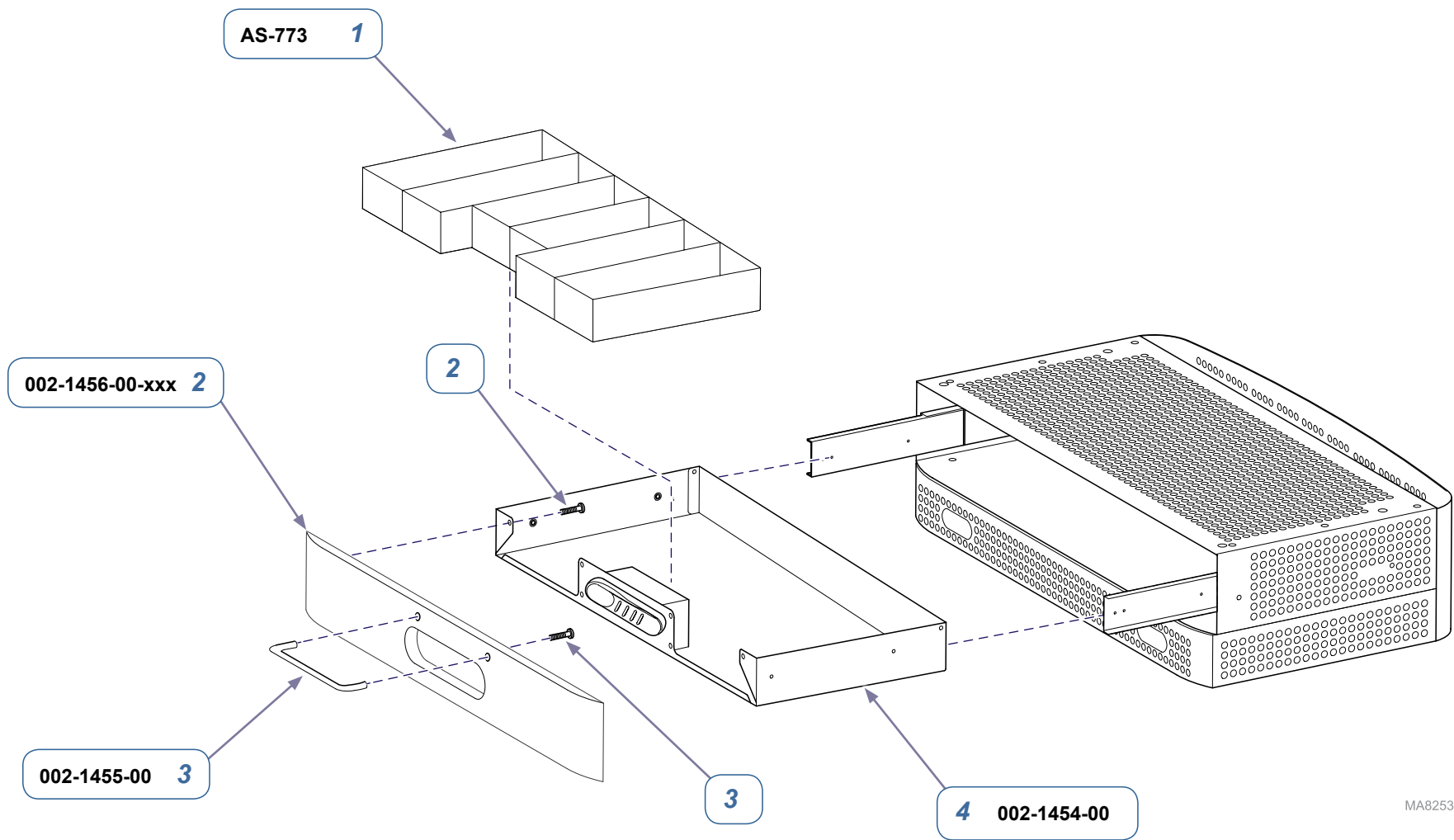

Models:	
Serial Numbers:	

Med Drawer Storage /w Lock		
Item	Description	Qty.
1	Med Drawer Insert	1
2	Kit, Drawer Front w/ Lock Cutout	1
3	Kit, Drawer Bin Pull	1
4	Kit, Locking Med Drawer Bin	1

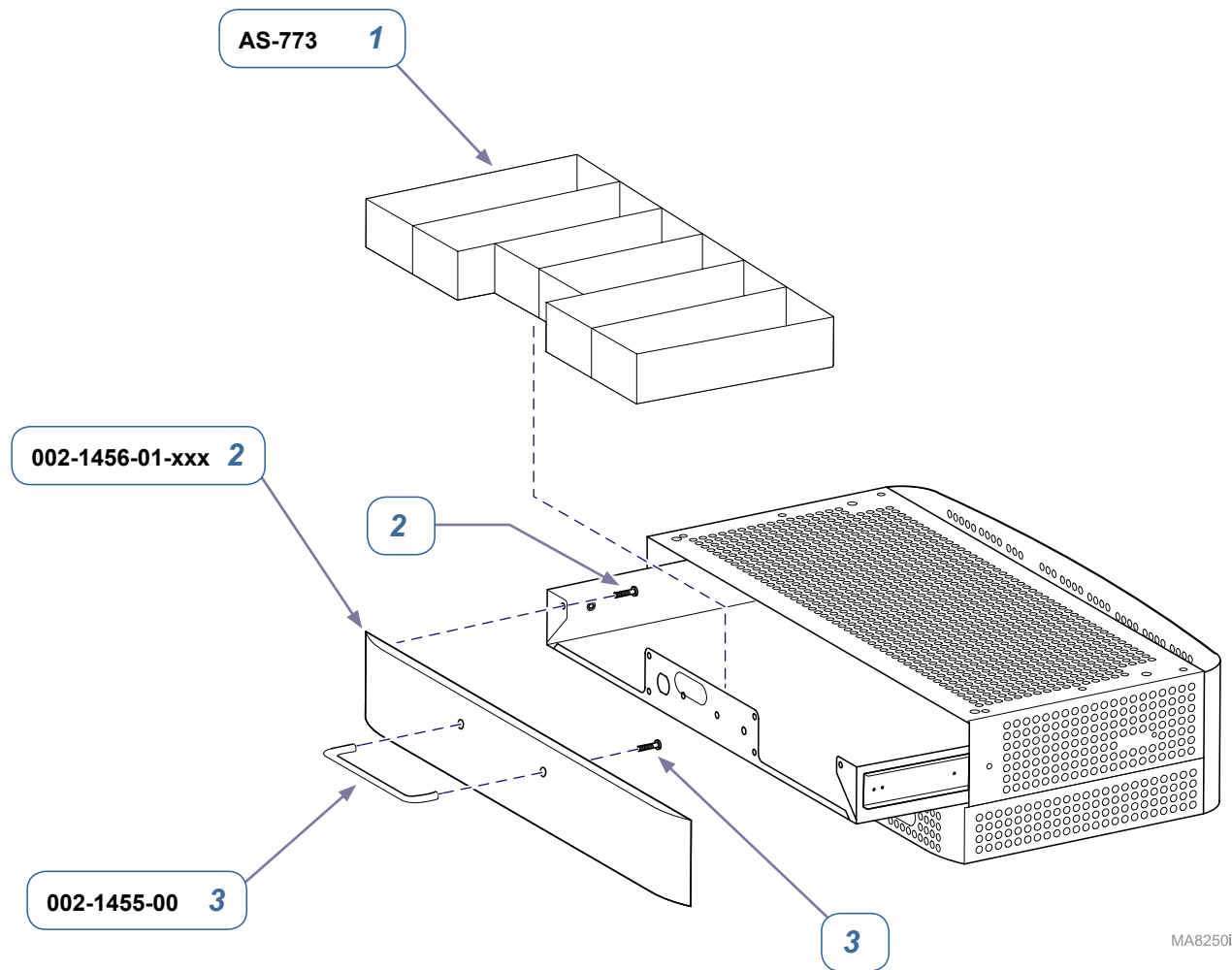

Models:	
Serial Numbers:	

Med Drawer Storage no Lock		
Item	Description	Qty.
1	Med Drawer Insert	1
2	Kit, Drawer Front w/o Lock Cutout	1
3	Kit, Drawer Bin Pull	1

SUBJECT TO CHANGE WITHOUT NOTICE

Midmark Corporation
For contact information, go to:
www.midmark.com

